

The evidence is in...

Enhanced Outcomes from Enhanced Stability

 **smith&nephew**
TRIGEN[®]
INTERTAN[®]
Intertrochanteric Antegrade Nail

Supporting healthcare professionals



Lower
risk of implant
failure and
non-union



Reduced
postoperative
pain



Faster
time to
fracture union



Proven
high return
to pre-fracture
status



How satisfied are you with current hip fracture outcomes?




Around **1 in 4** hip fracture patients over the age of 65 die within 12 months¹




Around **6.6%** will require reoperation due to complications²


And for those who survive:

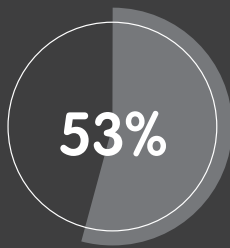


decline in mobility³

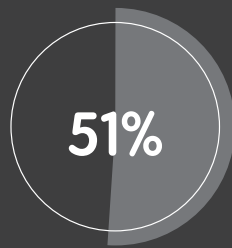


decline in daily living activities³

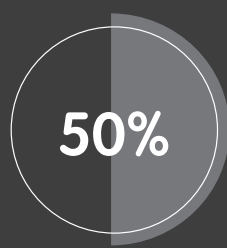

live with **reduced** ability to walk following fracture union¹

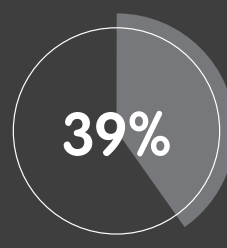

decline in self-reported health³

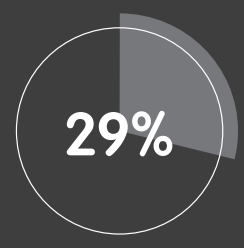

decline in fine motor skills³

Your patients can enjoy life after a hip fracture

The evidence is in! Based on data from more than two-dozen published studies, the TRIGEN[®] INTERTAN[®] Intertrochanteric Antegrade Nail allows patients to experience:



Lower risk of implant failure and non-union



Reduced postoperative pain




Faster time to fracture union



Proven high return to pre-fracture status

Here's how it works

"Success rate of the operation partly depends on factors that the surgeon cannot influence. Surgeons should therefore be aware of the factors that they can manipulate with a positive outcome." – Brujin et al, 2012

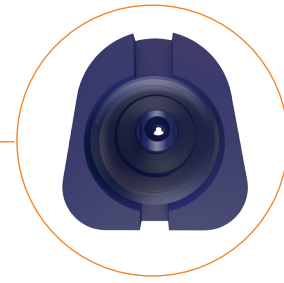


Intertrochanteric rotational stability

The trapezoidal shape provides a pressfit in the metaphyseal region and positions more material on the lateral side of the nail where tensile/stretching forces tend to be greatest




Maintain compression and eliminate Z-effect

Integrated Compression Screws thread together to generate push/pull forces that hold compression after instruments are removed and eliminate Z-effect

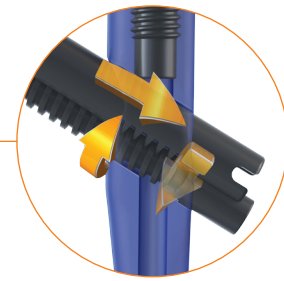

Control rotation during reduction

A worm gear mechanism converts rotation to active compression while stabilizing the medial fragment

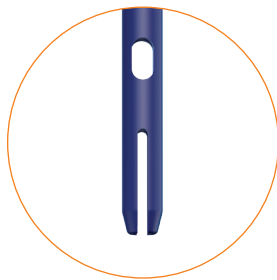

Prevent periprosthetic fractures

A clothes pin distal tip is less rigid to decrease the stress riser and reduce the incidence of anterior thigh pain



Eliminate medial migration

The head of the compression screw pushes medially against the nail and unloads stress forces off the lateral wall







Challenge Postoperative complications

Complication rates are still above 4% and can reach up to 16% in highly unstable fractures⁵



The TRIGEN[°]
INTERTAN[°] Solution:


**Lower risk of
implant failure and
non-union**^{6,7,8,9,10,11,12,13,14,15}

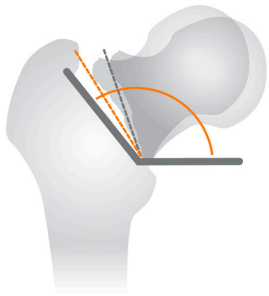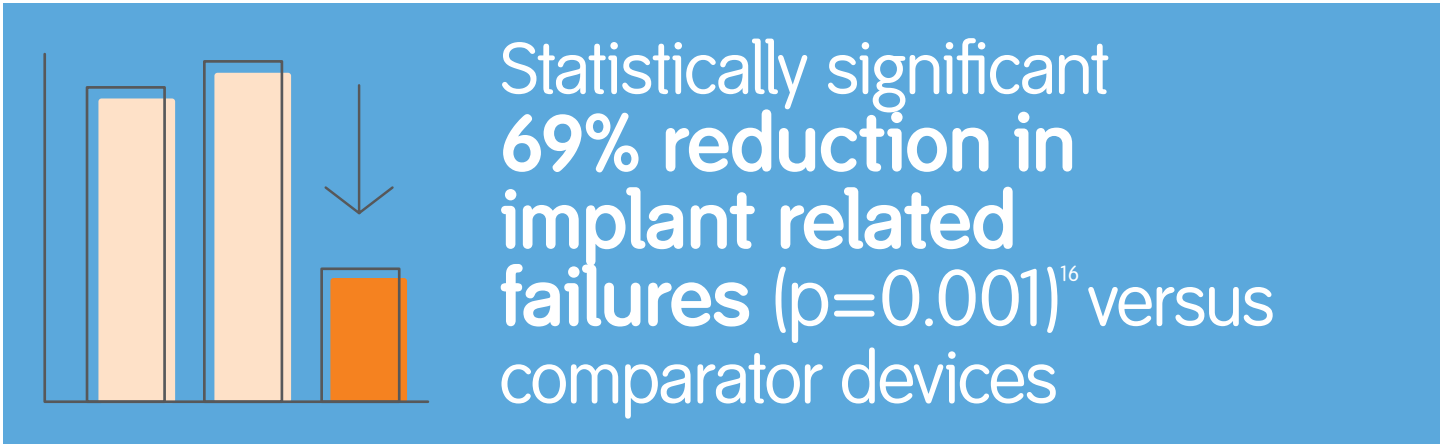

2.5x less varus

collapse when
compared to the
single screw¹⁷

- Less varus collapse^{5,17,18}
- Less peri-implant fractures^{13,19}
- Effective in reducing the potential role of the tip of the short nail as a stress riser¹³

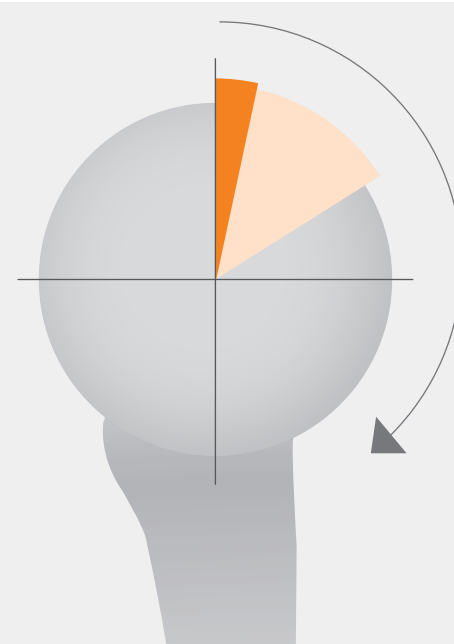

84% less


initial rotation
and varus
collapse⁵

Why INTERTAN?

The Integrated Compression Screws (ICS) of INTERTAN provide a second point of fixation in the femoral head, and allow for mechanical compression through the implant which is actively maintained after instrument removal. This combination creates strong interfragmentary friction and increases construct stability to resist complications such as rotation and varus collapse.



The TRIGEN[®] INTERTAN[®] Solution:
**Lower risk of implant failure
and non-union**^{6,7,8,9,10,11,12,13,14,15}

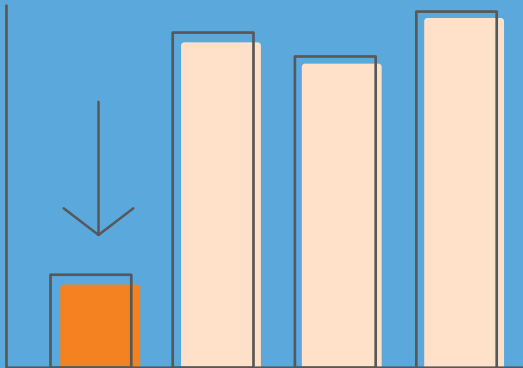

Statistically
significant **73%**
reduction in
non-union
($p=0.01$)¹⁶ versus
comparator devices

“The integrated dual screw device offered significantly increased stability throughout the time interval that would be needed for fracture healing.”

– Santoni et al, 2016



No non-unions in
radiographic analysis
of the TRIGEN
INTERTAN nail^{13,19,20}




Challenge pain management

Instability of the
bone-implant construct >
movement at the fracture
site > pain¹¹



The TRIGEN[®]
INTERTAN[®] Solution:
Reduced
postoperative
pain^{9,10,11,13,14,15}

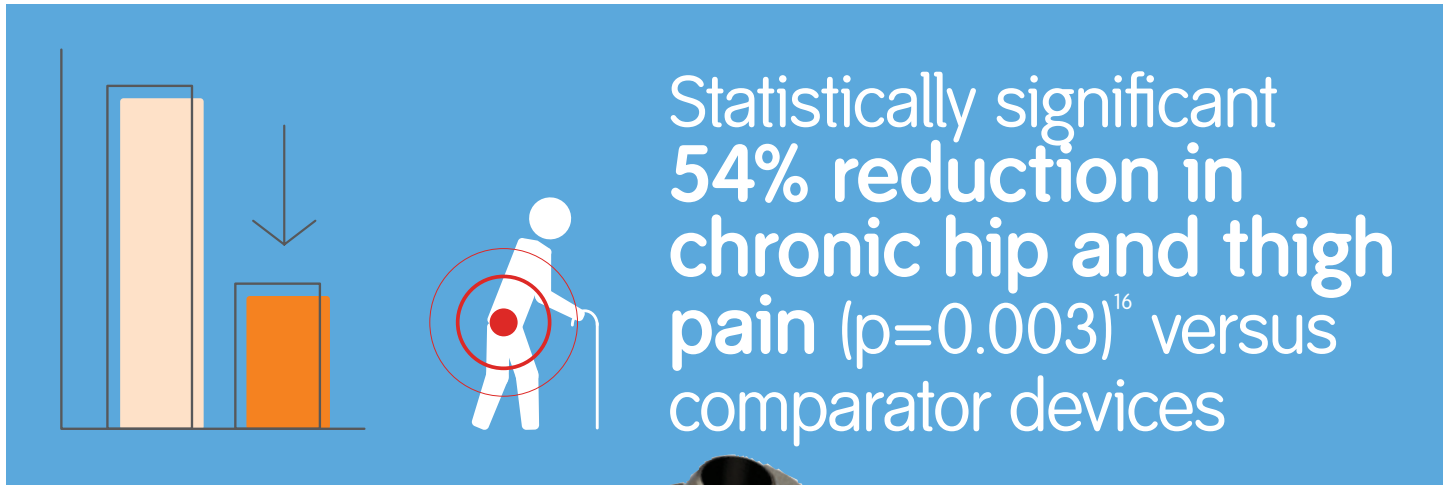

“When pain is not effectively managed, patients are not able to walk as they did before their injury, and they are more likely to have compromised pulmonary and cardiac function.” – Zanzone et al, 2016

“Poorly managed postoperative pain is associated with delayed ambulation, pulmonary complications, and delayed transition to lower levels of care.”
– Abou-Setta et al, 2011

Statistically significant
improvement in postoperative
pain and mobility¹⁹

“In our series, intertrochanteric fracture fixation using an INTERTAN nail lead to significantly shorter hospital stay, better functional outcomes, and less pain at 6 months.” – Berger-Groch et al, 2016



Why INTERTAN?

With compression actively maintained postoperatively using the ICS screws, INTERTAN is designed to reduce unnatural movement of the hip at the fracture site. Patients with INTERTAN have been shown to experience less pain and therefore may feel more comfortable weight bearing on their implant postoperatively.



Challenge Delayed healing

Insufficient stabilization > excessive motion of the fracture site > delayed healing



The TRIGEN[®] INTERTAN[®] Solution:

Faster time to fracture union^{9,11,12,13,19,}

23,24,25,26,27,28,29

Nearly 3 week faster time to fracture union¹⁶ versus comparator devices

TRIGEN INTERTAN

14.1 weeks

comparator

16.9 weeks




“Excess interfragmentary shear or rotational movements inhibit repair and can result in a significant delay to healing.” – Gaston et al, 2007

“The INTERTAN nail can reduce healing time and is a good choice for elderly patients who need to walk bearing full weight in the early post-operative period.” – Zhang et al, 2013

Simulated gait

Simulated chair rise

Simulated chair rise

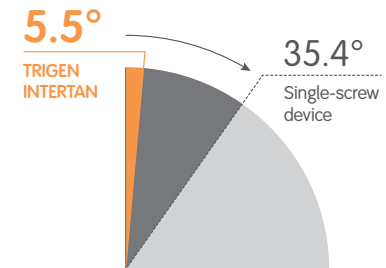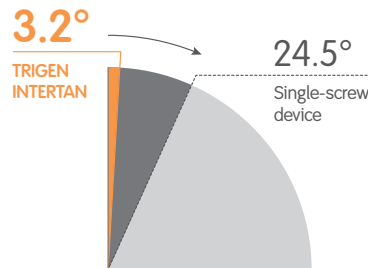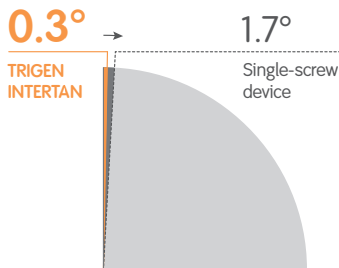

5x greater

initial rotational stability⁵

In a biomechanical simulated gait study comparing TRIGEN INTERTAN and Gamma3

7x less

femoral head rotation¹⁷

In a biomechanical simulated chair rise study comparing TRIGEN INTERTAN and Gamma3

7x reduction

in maximum femoral head rotation¹⁷

In a biomechanical simulated chair rise study at the end of 4x body weight loading or until failure

Why INTERTAN?

By properly stabilizing the anatomy and maintaining an anatomical reduction, The INTERTAN ICS screws resist excessive motion in order to create a more stable healing environment. This provides the patient’s biology a better chance to achieve an earlier and more successful union at the fracture site.



Challenge Poor functional outcomes

Femoral neck shortening
> decreases moment arm
of abductors > reduced
patient function³¹




The TRIGEN[®]
INTERTAN[®] Solution:
**Proven high
return to
pre-fracture
status**^{6,13,19,23}

Statistically significant **higher SF-36 score** in favor of TRIGEN INTERTAN ($p=0.002$)¹⁶ versus the comparator in one study⁶



up to 93% return to pre-fracture ambulatory status

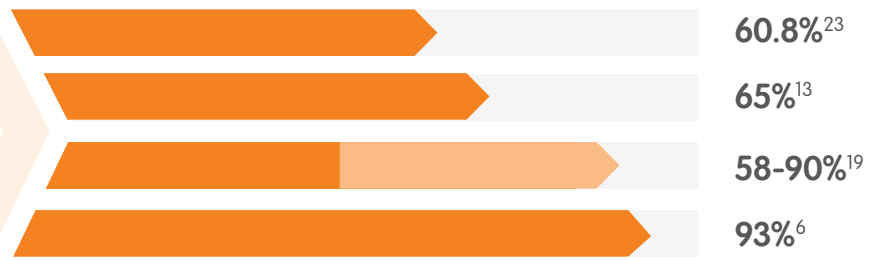

“A large proportion of those patients who survive never recover to their prefracture level of function.”
– Abou-Setta et al, 2011

“Shortening of greater than 2cm is known to adversely affect locomotor function in otherwise active individuals.” – Sanders et al, 2017

“Shortening of the femoral neck was the only significant variable predictive of a low SF-36 physical functioning score.” – Zlowodzki et al, 2008



No uncontrolled collapse of the neck¹⁹

Less femoral neck shortening^{8,9,17}


Why INTERTAN?

Utilizing the ICS screws for controlled active compression, rather than relying on weight bearing and uncontrolled sliding, helps resist shortening of the femoral neck which can improve patient function. By restoring the patient’s natural anatomic measurements and preserving limb length, INTERTAN results in highly successful postoperative ambulatory outcomes.


Inadequate
compression

Single Screw

“Because the screw is rotationally unstable within the bone when using a single lag screw, flexion-extension of the limb results in loosening of the bone-screw interface, with the screw secondarily cutting out” – Zhang et al, 2013

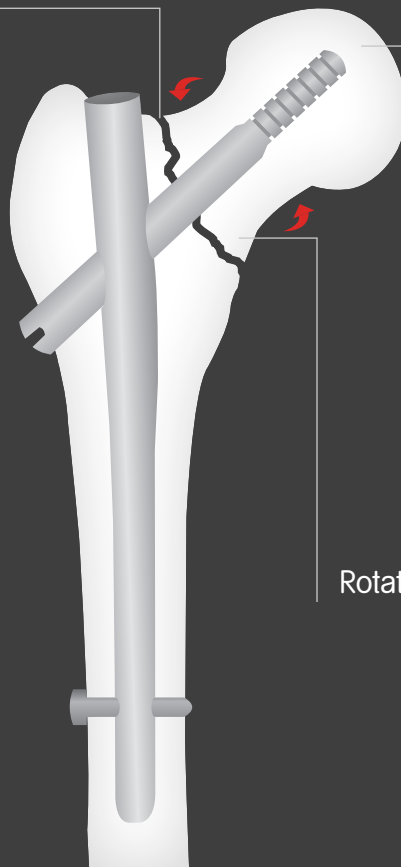

Single Pivot Point

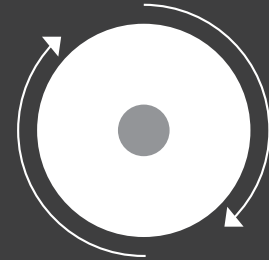


Rotational Instability

Compression
Maintained

Integrated Compression Screws

“With the more recent identification that rotational instability contributes to malunion and implant-bone construct failure, the use of an integrated-slide implant should be considered to provide added rotational stability in unstable fracture patterns.” – Baldwin et al, 2016

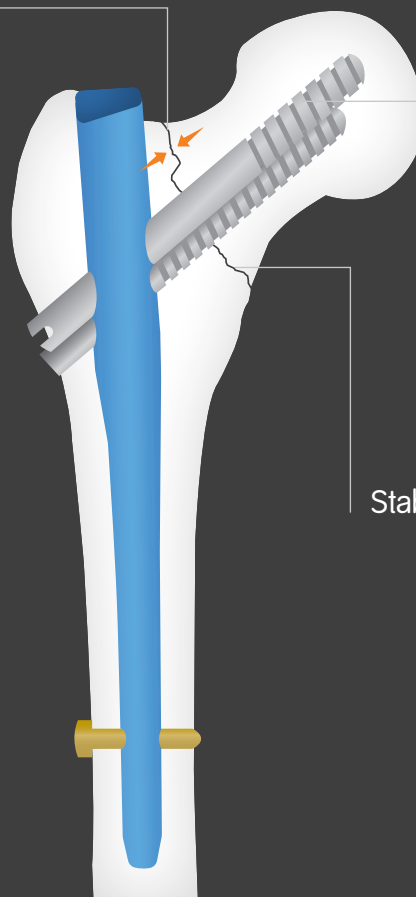


Two Points of Fixation



Stable Fracture Site

There is a lack of definitive evidence identifying the superiority of a helical blade or lag screw implant³³

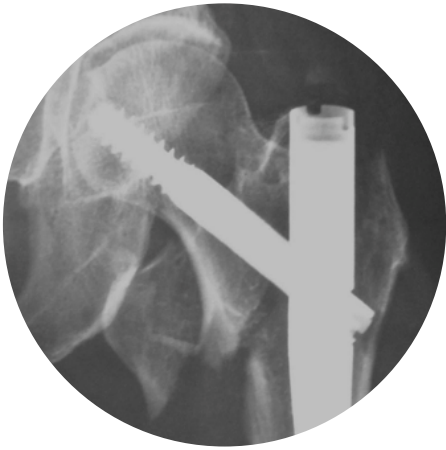

“The use of the INTERTAN system may be an improvement in surgery compared to Gamma 3...In our department, we have standardized the use of INTERTAN nail for the treatment of intertrochanteric fractures.” - Su et al, 2016

“The results of our study show that the incidence of femoral shaft fractures, rotational loss of reduction, varus collapse of the head/neck, [...] cut-out, and femoral neck shortening were decreased in group IT comparing with group PFNA-II.”
- Yu et al, 2016

“INTERTAN is superior to DHS in internal fixation stability, thus better applies in cases of osteoporosis and unstable fractures.”
- Wang et al, 2014



Economic Impact:

In the changing economic landscape of healthcare, better patient outcomes mean better outcomes for hospitals. When considering the costs involved in treating a patient with a hip fracture, the benefits of the TRIGEN[®] INTERTAN[®] system – **lower risk of implant failure and non-union, reduced postoperative pain, faster time to fracture union, and a proven high return to pre-fracture status** – can help you achieve better outcomes more efficiently.



“The priority remains improving functional outcomes and reducing complications. If, as a profession, we are to rise to the challenge of the ageing population, more is going to be needed for less.”

–Ollivere et al, 2017



Smith & Nephew, Inc.

1450 Brooks Road
Memphis, TN 38116
USA

Telephone: 1-901-396-2121
Information: 1-800-821-5700
Orders/Inquiries: 1-800-238-7538

www.smith-nephew.com

®Trademark of Smith & Nephew.
©2017 Smith & Nephew.
05036 V2 0817

Supporting healthcare professionals for over 150 years

References

1. American Academy of Orthopaedic Surgeons. Hip fractures in seniors: a call for health system reform. Position Statement 1144. Rosemont, IL: 1999. **2.** Mundi S et al. Similar mortality rates in hip fracture patients over the past 31 years: A systematic review of RCTs. *Acta Orthopaedica* 2014; 85(1): 54-59. **3.** Bentler SE, Liu L, Obrizan M, Cook EA, Wright KB, Geweke JF, et al. The aftermath of hip fracture: discharge placement, functional status change, and mortality. *Am. J Epidemiol.* 2009 Nov 15;170(10):1290-9. doi: 10.1093/aje/kwp266. **4.** Brujin K, Hartog D, Tuinebreijer W, Roukema G. Reliability of Predictors for Screw Cutout in Intertrochanteric Hip Fractures. *J Bone Joint Surg Am.* 2012;94:1266-1272. **5.** Hoffmann S, Paetzold R, Stephan D, Püschel K, Buehren V, Augat P. Biomechanical evaluation of interlocking lag screw design in intramedullary nailing of unstable peritrochanteric fractures. *J Orthop Trauma.* 2013;27(9):483-490. **6.** Berger-Groch J, Rupperecht M, Schoepper S, Schroeder M, Rueger JM, Hoffmann M. Five-Year Outcome Analysis of Intertrochanteric Femur Fractures: A Prospective Randomized Trial Comparing a 2-Screw and a Single-Screw Cephalomedullary Nail. *J Orthop Trauma.* 2016;30:483-488. **7.** Matre K, Vinje T, Havelin LI, et al. TRIGEN INTERTAN intramedullary nail versus sliding hip screw: a prospective, randomized multicenter study on pain, function, and complications in 684 patients with an intertrochanteric or subtrochanteric fracture and one year of follow-up. *J Bone Joint Surg Am.* 2013;95:200-208. **8.** Sanders D, Bryant D, Tieszer C, et al. A Multicenter Randomized Control Trial Comparing a Novel Intramedullary Device (InterTAN) Versus Conventional Treatment (Sliding Hip Screw) of Geriatric Hip Fractures. *J Orthop Trauma.* 2017;31:1-8. **9.** Seyhan M, Turkmen I, Unay K, Ozkut AT. Do PFNA devices and Intertan nails both have the same effects in the treatment of trochanteric fractures? A prospective clinical study. *J Orthop Sci.* 2015;20:1053-1061. **10.** Su H, Sun K, Wang X. A randomized prospective comparison of Intertan and Gamma3 for treating unstable intertrochanteric fractures. *Int J Clin Exp Med.* 2016;9:8640-8647. **11.** Zhang S, Zhang K, Jia Y, Yu B, Feng W. InterTan nail versus Proximal Femoral Nail Antirotation-Asia in the treatment of unstable trochanteric fractures. *Orthopedics.* 2013;36:e288-294. **12.** Wang Q, Yang X, He HZ, Dong LJ, Huang DG. Comparative study of InterTAN and Dynamic Hip Screw in treatment of femoral intertrochanteric injury and wound. *Int J Clin Exp Med.* 2014;7:5578-5582. **13.** Wu Y, Watson JT, Kuldjanov D, Jackman J. Rotationally stable fixation for intertrochanteric hip fractures: the Intertan experience, surgical technique, and outcomes. *Techniques in Ortho.* 2014;29:3:120-132. **14.** Yu W, Zhang X, Zhu X, Hu J, Liu Y. A retrospective analysis of the InterTan nail and proximal femoral nail anti-rotation-Asia in the treatment of unstable intertrochanteric femur fractures in the elderly. *J Orthop Surg Res.* 2016;11:10. **15.** Zehir S, Sahin E, Zehir R. Comparison of clinical outcomes with three different intramedullary nailing devices in the treatment of unstable trochanteric fractures. *Ulus Travma Acil Cerrahi Derg.* 2015;21(6):469-476. **16.** Leo N, Dunbar C, Ridgway J, Horner A. The TRIGEN INTERTAN Intertrochanteric Antegrade Nail: A Systematic Literature Review and Meta-analysis of Clinical Outcomes Compared to Standard of Care in the Treatment of Intertrochanteric Hip Fractures. *Bone&Joint Outcome.* 2017;4(1):1-20. Lit no: 10205 V1 07/17. **17.** Serrano-Riera R, Blair JA, Downes K, Sanders R. Cephalo-medullary nail fixation of intertrochanteric fractures: are two proximal screws better than one? Abstract presented at: Orthopaedic Trauma Association Annual Meeting; October 15-18, 2014; Tampa, FL, USA. **18.** Santoni B, Nayak A, Cooper S, et al. Comparison of Femoral Head Rotation and Varus Collapse Between a Single Lag Screw and Integrated Dual Screw Intertrochanteric Hip Fracture Fixation Device Using a Cadaveric Hemi-Pelvis Biomechanical Model. *J Orthop Trauma.* 2016;30:164-169. **19.** Ruecker AH, Rupperecht M, Gruber M, Gebauer M, Barve s using an intramedullary nail with integrated cephalocervical screws and linear compression. *J Orthop Trauma* 2009;23:22-30. **20.** Galli M, Ciriello V, Bocchino L, Gangemi NM, Peruzzi M, Marzetti E. Clinical and functional outcomes of internal fixation with intertrochanteric antegrade nail in older patients with proximal extracapsular femoral fractures. *Eur J Trauma Emerg Surg.* 2013;10/17 2013:1-6. **21.** Zanzone A. Current Challenges in Pain Management in Hip Fracture Patients. *J Orthop Trauma.* 2016;30:S1-S5. **22.** Abou-Setta A, Beaupre L, Jones C, et al. Pain Management Interventions for Hip Fracture. Agency for Healthcare Research and Quality. 2011;Publication No. 11-EHC022-EF. **23.** Kim JW, Kim TY, Ha YC, Lee YK, Koo KH. Outcome of intertrochanteric fractures treated by intramedullary nail with two integrated lag screws: A study in Asian population. *Indian J Orthop.* 2015;49:436-441. **24.** Tao R, Lu Y, Xu H, Zhou ZY, Wang YH, Liu F. Internal fixation of intertrochanteric hip fractures: a clinical comparison of two implant designs. *Sci World J.* 2013;2013:1-8. **25.** Huang FT, Lin KC, Yang SW, Renn JH. Comparative study of the proximal femoral nail antirotation versus the reconstruction nail in the treatment of comminuted proximal femoral fracture. *Orthopedics.* 2012;35:e41-47. **26.** Sahin EK, Imerci A, Kinik H, Karapinar L, Canbek U, Savran A. Comparison of proximal femoral nail antirotation (PFNA) with AO dynamic condylar screws (DCS) for the treatment for unstable peritrochanteric femoral fractures. *Eur J Orthop Surg Traumatol.* 2014;24:347-352. **27.** Hsueh K, Fang C. Risk factors in cutout of sliding hip screw in intertrochanteric fractures: an evaluation of 937 patients. *Int Orthop.* 2010;34:1273-1276. **28.** Liu Y, Tao R, Liu F, et al. Mid-term outcomes after intramedullary fixation of peritrochanteric femoral fractures using the new proximal femoral nail antirotation (PFNA). *Injury.* 2010;41:810-817. **29.** Riha D, Bartoniček J. Internal fixation of peritrochanteric fractures using DHS with a two-hole side-plate. *Int Orthop.* 2010;34:877-882. **30.** Gaston MS, Simpson AHRW. Inhibition of fracture healing. *J Bone Joint Surg Br.* 2007;89-B:1553-1560. **31.** Rueger J, Moore C. Shortening of the femoral neck following peritrochanteric fracture. *Bone Joint Sci.* 2011 May;2(5). **32.** Zlowodzki M, Brink O, Switzer J, et al. The effect of shortening and varus collapse of the femoral neck on function after fixation of intracapsular fracture of the hip. *J Bone Joint Surg Br.* 2008;90:1487-1494. **33.** Baldwin P, Lavender R, Sanders R, Koval K. Controversies in Intramedullary Fixation for Intertrochanteric Hip Fractures. *J Orthop Trauma.* 2016;30:635-641. **34.** Öllivere B, Das A, Shiyji F. Hip fractures: The state of the art in 2017. *The British Editorial Society of Bone & Joint Surgery.* 2017;6:3.