



- Shoulder pleats allow freedom of movement for the upper body.

- Forward facing sleeves correspond with natural working position to allow comfort and freedom of movement.


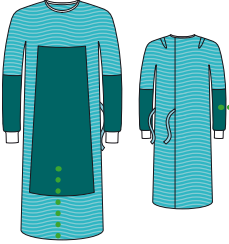


- Long sleeve increases freedom of movement and reduces risk of glove sliding.

BARRIER® Surgical Gown - ULTIMATE

New

- **Highly fluid repellent** soft SMS material for advanced protection
- **Treated** to provide increased repellency to low surface tension fluids
- **Very soft** and breathable for ultimate comfort
- **Ergonomic design** for comfort and freedom of movement
- **Suitable** for all procedures.

BARRIER Surgical Gown - ULTIMATE

	Ref. No.	Description	Size	Disp/Trp
Standard Performance	STANDARD PERFORMANCE			
	680102	SG SP two towels, wrapped	L 121 cm	18/36
	680103	SG SP two towels, wrapped	L-L 137 cm	16/32
	680104	SG SP two towels, wrapped	XL 132 cm	16/32
	680105	SG SP two towels, wrapped	XL-L 148 cm	16/32
	680106	SG SP two towels, wrapped	XXL-L 147 cm	14/28
	HIGH PERFORMANCE			
  	690102	SG HP two towels, wrapped	L 121 cm	16/32
	690103	SG HP two towels, wrapped	L-L 137 cm	14/28
	690104	SG HP two towels, wrapped	XL 132 cm	14/28
	690105	SG HP two towels, wrapped	XL-L 148 cm	14/28
	690106	SG HP two towels, wrapped	XXL-L 147 cm	12/24
	690107	SG HP two towels, wrapped	XXL-XL 167 cm	12/24
	SG = Surgical Gown SP = Standard Performance HP = High Performance			

- Impermeable front reinforcement to minimise risk of strike through
- Breathable for comfort.
- Impermeable sleeve reinforcement extended well above elbow to minimise risk of strike through
- Ultrasonic welded sleeve for increased repellency
- Absorbing and soft for comfortable feeling.

Gown materials	Gown Body & Sleeve Front Reinforcement Sleeve Reinforcement	Nonwoven pulp/polyester fibres Polyethylene plastic film (breathable) Polyethylene plastic film, viscose/polyester nonwoven
Product classification	EN 13795 Standard Performance and AAMI Level 2 & EN 13795 High Performance and AAMI Level 3	
Packaging	Single packed gowns provided in a double box packaging system, ensures protection against dust, moisture and damage during transport, handling and storage.	
Environment	At Mölnlycke Health Care we have a long-standing commitment to minimise the impact of our activities on the environment, without compromising the efficiency and reliability of our products and services. Sound environmental practice is a priority in our production and at our sites around the world.	

Find out more at www.molnlycke.co.uk

Mölnlycke Health Care Ltd, The Arenson Centre, Arenson Way, Dunstable, Bedfordshire, LU5 5UL.
T 0870 6060766. www.molnlycke.co.uk

The Mölnlycke Health Care name and logo and BARRIER® and ProcedurePak® are registered trademarks of Mölnlycke Health Care AB. Copyright (2012)