


MEDICAL DEVICES

Injectomat TIVA Agilia

ANAESTHESIA SYRINGE PUMP

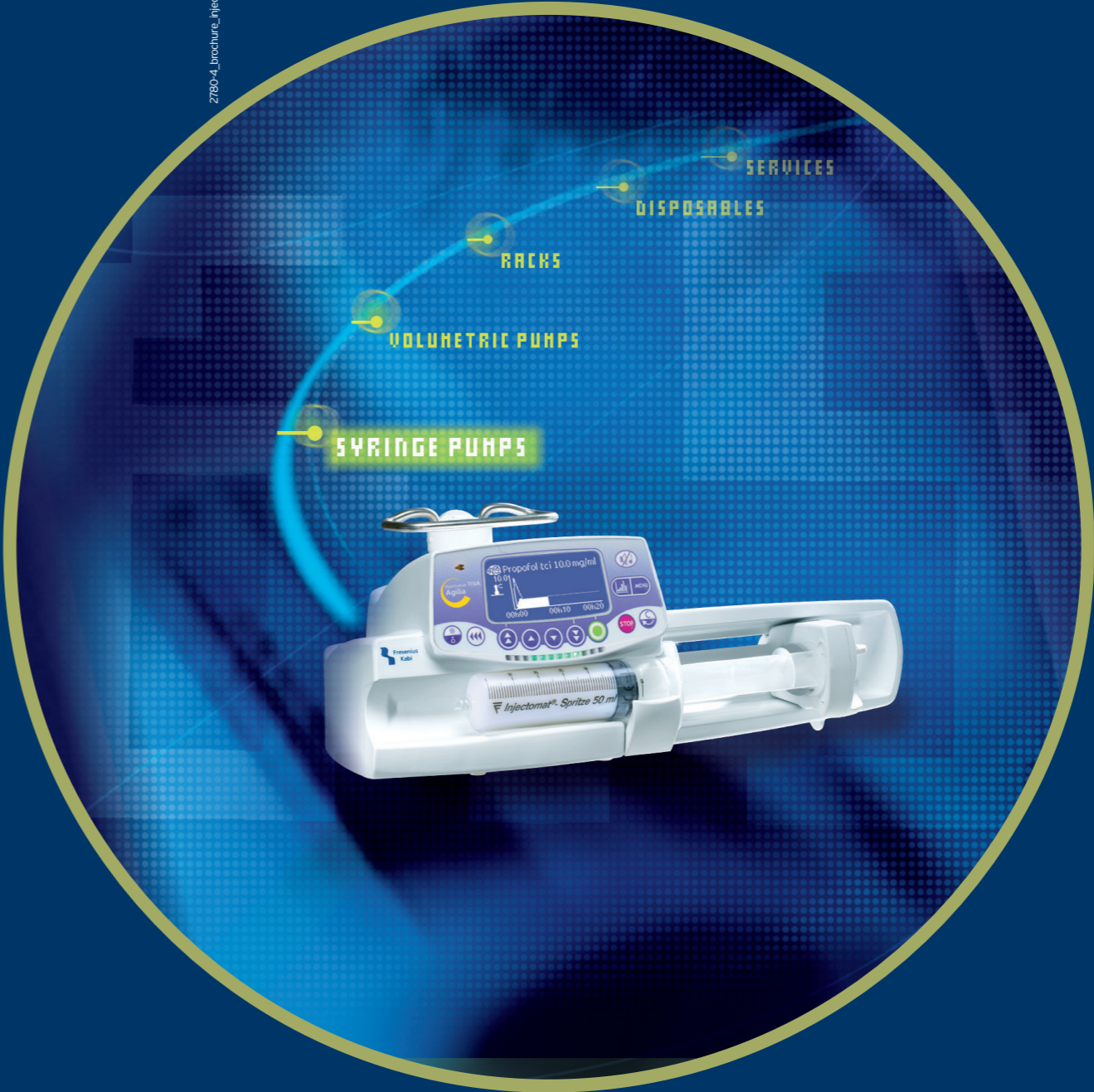

Z780-4\_brochure\_injectomat\_tiva\_agilia\_eng

# Intuitive intravenous anaesthesia

Injectomat TIVA Agilia is a single channel syringe pump dedicated to anaesthesia.

Depending on your medical practice, it can be used either in TIVA\* mode (ml/h or dose rate), or in TCI\* mode (plasma concentration or effect site concentration).

The device combines ease of use and a complete set of features, including PK\* models for 4 drugs: propofol, remifentanyl, sufentanil, alfentanil (in option).

Like all the Agilia devices, Injectomat TIVA Agilia is highly intuitive to use. The learning process is simple and short thanks to a clear front panel, obvious keys and meaningful pictograms.

Numerous safety features guarantee a secured use of the device: syringe positioning detector, clear texts and diagrams, and advanced pressure management. Combined with Vigilant Drug'Lib, intravenous anaesthesia becomes even safer.

\* TIVA: Total Intravenous Anaesthesia


\* TCI: Target Controlled Infusion

\* PK: Pharmacokinetic




## Intuitive confidence

The infusion setting is guided and each step clearly identified. The displayed instructions are simple to follow. Intuitive settings increase user confidence.

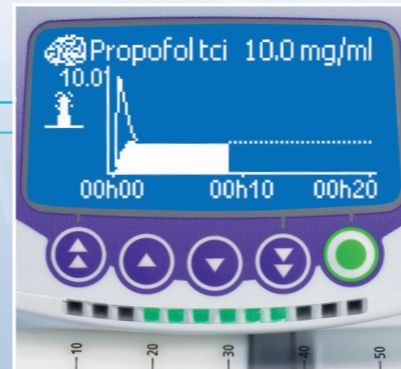

## Obvious monitoring

Monitoring is simple thanks to easily accessible numerical and graphical information. The syringe is always visible and the flow indicator lights enable a quick and safe monitoring.

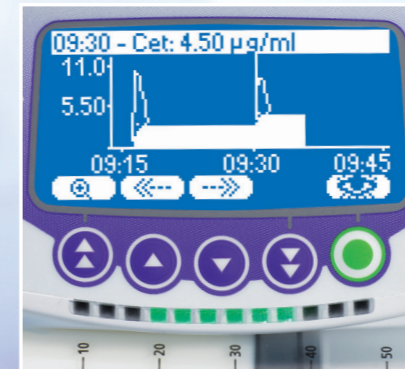

## Quick history access



The graph key gives direct access to the:

- Pressure history,
- Dose rate history,
- Target and concentration history,
- Infused dose history.



## Vigilant, the IV Medication Safety Solution

Assisted with Vigilant Drug'Lib, Injectomat TIVA Agilia brings you a higher level of safety in IV drug administration.

Vigilant Drug'Lib is a software programme that enables the creation of personalised drug protocols. Each drug is defined with a default dilution, a prescription unit and a range of authorised settings.

In the Vigilant Drug'Lib mode the infusion setting is simplified and shortened. The risks of under or overdosage of IV drugs are limited.

## Agilia

The Agilia family offers a complete range of infusion pumps, accessories and disposables to cover all major IV drug administration needs.

### Injectomat Agilia



### Injectomat MC Agilia



### Injectomat TIVA Agilia



### Volumat Agilia



### Volumat MC Agilia



### Links & Accessories



### Administration Sets



### Maintenance Tools



Injectomat, Volumat, Agilia and Vigilant are registered trademarks by Fresenius Kabi in selected countries.