

# CODAN Chemoprotect®



Complete protection in  
preparation of cytotoxic drugs




# CODAN Chemoprotect®



The CODAN Chemoprotect® product range offers safe and reliable protection of personnel in the handling of cytotoxic drugs. The combined use of specially matched individual products considerably minimizes the risk of contaminations in the regular handling of cytotoxic drugs, as well as in the event of accidents (spills). Chemoprotect® is the perfect addition to CODAN CYTO® products for the preparation and administration of cytotoxic drugs.

The exceptional reliability of the Chemoprotect® system is confirmed by the certification of the most important products as Personal Protective Equipment (PPE) according to the European Directive 89/686/EEC.

## Chemoprotect® Latex-free protective gloves

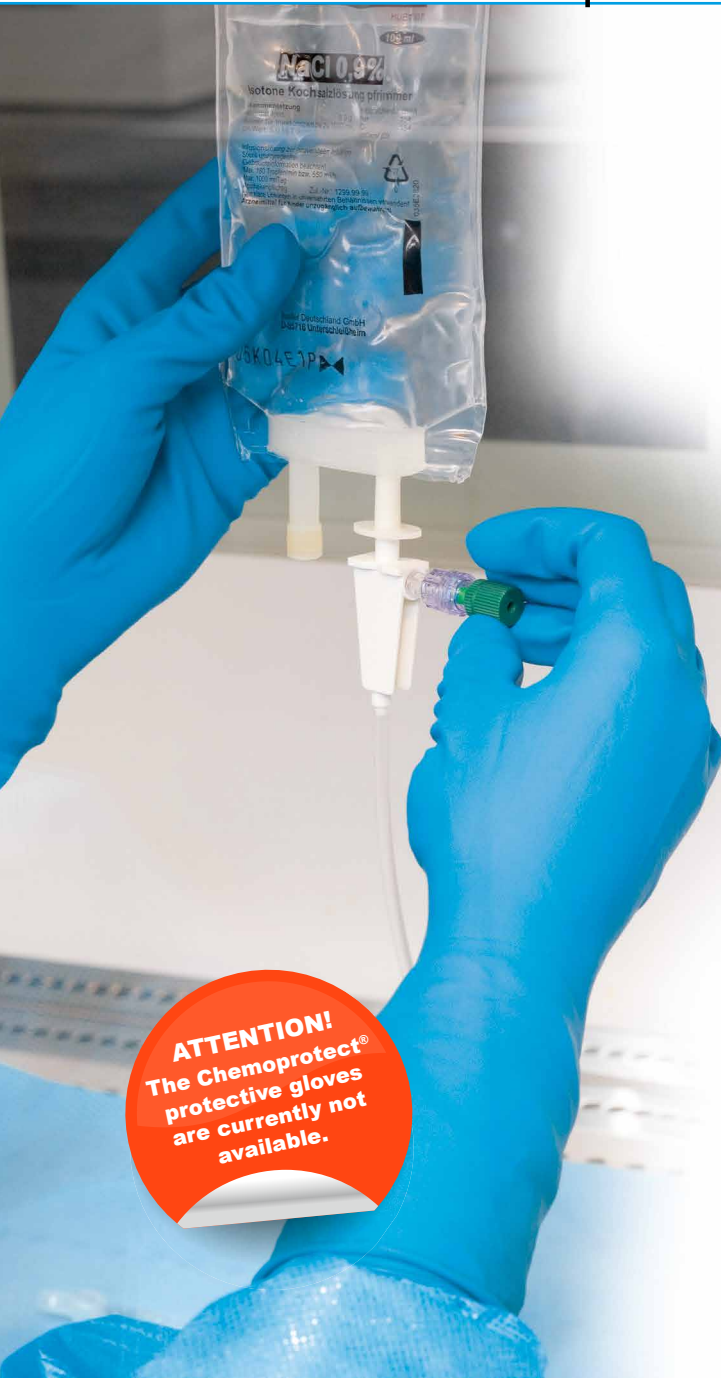

### Precise work with maximum protection

The new latex-free Chemoprotect® protective gloves offer reliable protection in the handling of cytotoxic drugs. They are distinguished by their close fit and excellent sensory properties resulting from different thicknesses of the thin material in the various areas of the gloves.

These Latex-free gloves are made of neoprene (polychloroprene) and are provided with an optimal chemical barrier. As a result, they also offer allergy sufferers excellent wearing convenience. Despite their long cuffs, a special inner lining allows easy fitting and removal of the gloves even without powder.

Chemoprotect® protective gloves are certified as Personal Protective Equipment (PPE) according to the European Directive 89/686/EEC. They are available sterile or unsterile in four glove sizes from XS to L.

- No allergy risk
- Optimal chemical barrier
- Excellent sensory properties
- Well-fitting
- Long glove cuffs
- Fast and easy fitting and removal
- Polycoated inner lining
- Individually packed in pairs either sterile or unsterile
- Dispenser appropriate for wards

**ATTENTION!**  
The Chemoprotect®  
protective gloves  
are currently not  
available.



# Chemoprotect® Prep mat

## Safe base

The Chemoprotect® prep mat has been specially developed for cytotoxic drug preparation. The three-layer structure of the material allows safe working and risk-free subsequent disposal. The prep mat is available sterile or non-sterile in two sizes for most workplace dimensions.

Dimensions:

Large prep mat: 41 x 54 cm

Small prep mat: 27 x 42 cm



STERILE EO

Lint-free  
fluid-permeable surface

Highly absorbent  
fluid-retaining middle layer

Fluid-impermeable  
bottom layer

# Chemoprotect® Arm sleeves

## Proven alternative

The preferred use of Chemoprotect® arm sleeves has proven successful in some workplaces. With regard to their material and processing properties they are identical to the Chemoprotect® protective gown. An elastic rubber cuff at the upper end provides a safe and comfortable fit with very good freedom of movement.



STERILE EO

CE 0299



The Chemoprotect® bag with a safe-seal closure is ideal for the transport of medication to the ward and is also excellent for the disposal of substances. The bags are available in transparent as well as in white (for light-sensitive cytotoxic drugs).



# Chemoprotect® Spill Box

## The complete accident set

With the Chemoprotect® Spill Box, CODAN offers a complete system for immediate and effective contamination protection in the event of accidents with cytotoxic drugs.

Alongside a complete Personal Protective Equipment (PPE) for the user, the Chemoprotect® Spill Box contains all necessary equipment for the disposal of contaminated objects and for the safe initial cleaning of contaminated surfaces.

The PPE products in the Chemoprotect® Spill Box consist of a Chemoprotect® protective gown, TouchNTuff® protective gloves in two sizes, sturdy over-gloves and overshoes, protective goggles along with a respiratory protection mask for reliable protection.

Barricade tape and warning signs are used to safeguard the accident location. A special powder binds the spilled substances which can then be removed with special cloths, wipers and shovels and can subsequently be disposed of together with all contaminated objects in specially labelled bags and containers. In the event of ampoule glass breakage, the debris and fluid can be removed with cotton wool and tweezers.

## Chemoprotect® Spill Box:

- Instructions for use
- Manufacturer's information
- Chemoprotect® protective gown
- TouchNTuff® protective gloves (2 sizes)
- Breathing mask
- Protective goggles
- Protective latex gloves (405)
- Overshoes (shoe covers)
- Emergency signs
- Chemoprotect® waste bag
- Powersorb cloths
- Disposal Box
- Wire binder
- Wiper
- Green-Z® powder
- Tongs
- Cotton wool
- Scoop
- Tape





# Chemoprotect® Protective gown

## Practical safety precautions

The greatest possible user protection in the handling of cytotoxic drugs is guaranteed by the Chemoprotect® protective gown. It is certified as Personal Protective Equipment (PPE) according to the European Directive 89/686/EEC.

The fluid-impermeable coating of the front of the gown and sleeves plus the extra heat-sealed seams provide reliable safety. The elasticated arm cuffs are hypoallergic (diazo-free) and offer a safe transition to the protective gloves.

The Chemoprotect® protective gown is available either sterile or non-sterile and in both cases in three sizes (S/M, L, XL). The lint-free material with good breathing properties and a comfortable cut ensure high wearing convenience.



| **STERILE EO**

**CE 0299**



# CODAN Worldwide



CODAN is known internationally as a manufacturer and supplier of disposable medical transfer systems. The CODAN Companies have more than 1500 employees around the world.

The name CODAN is synonymous with reliability, quality and precision based on the know-how and experience gained from more than 60 years of research and development. Company-owned production facilities and sales companies around the world are a guarantee for efficient production, a tight-knit sales network and a first-class service for customers in the healthcare sector.

## CODAN Companies

CODAN Medizinische Geräte GmbH & Co KG · Deutschland  
CODAN pvb Critical Care GmbH · Deutschland  
CODAN pvb Medical GmbH · Deutschland  
CODAN PORTUGAL, S.A. · Portugal  
CODAN 11, S.A. · Portugal  
CODAN US Corporation · California · USA  
CODAN Inc. · California · USA  
CODAN NORGE AS · Norge  
CODAN TRIPLUS AB · Sverige  
CODAN Limited · Great Britain  
CODAN FRANCE Sarl · France  
CODAN Medical AG · Schweiz  
CODAN ARGUS AG · Schweiz  
CODAN BV · Nederland  
CODAN s.r.l. · Italia  
CODAN Medical GmbH · Österreich  
CODAN Steritex ApS · Danmark  
CODAN Medical ApS · Danmark  
CODAN DEHA ApS · Danmark

### Manufacturer

CODAN Medizinische Geräte GmbH & Co KG  
Stig Husted-Andersen Straße 11 · 23738 Lensahn, Germany  
Tel. +49 (0) 43 63 · 51 11 · Fax. +49 (0) 43 63 · 51 12 14  
[www.codan.de](http://www.codan.de)

**The decisive connection**



### CODAN Product range

- Infusion sets
- Transfusion sets
- Extension lines and manifold connectors
- Infusion and transfusion accessories
- Infusion filters and filter systems
- Neonatology/Paediatric products
- Withdrawal, preparation and administration systems
- CODAN CYTO®
- Chemoprotect® products
- CODAN ENTERAL
- Single use syringes
- Invasive blood pressure monitoring systems
- Infusion pumps
- Continence care products
- Other CODAN Products

Compliance of the established quality management systems with the provisions of EN ISO 13485 and the Council Directive 93/42/EEC has been certified by the relevant, competent notified bodies. Our products are thus entitled to bear the following marks:

**CE 0123 TÜV SÜD Product Service GmbH**

CODAN Medizinische Geräte GmbH & Co KG  
23738 Lensahn, Germany

CODAN pvb Critical Care GmbH  
85661 Forstinning, Germany

CODAN US Corporation · Santa Ana, CA 92704, USA  
CODAN ARGUS AG · 6340 Baar, Switzerland

**CE 0543 Presafe Denmark A/S**

CODAN Steritex ApS · 3060 Espergærde, Denmark  
CODAN Medical ApS · 4970 Rødby, Denmark

