

clinell[®]

www.clinell.com

**NEW
FORMULA
STRONGER
CLEANING
POWER, ULTRA
LOW SMEAR**

CUTTING EDGE TECHNOLOGY COMBINED WITH EASE OF USE TO CREATE AN INCREDIBLY EFFECTIVE AND ENVIRONMENTALLY FRIENDLY FORMULA

Unique new patent pending formula

Natural and fully biodegradable, with ultra low smear.

Optimum cleaning

Even stronger cleaning power, developed specifically for use in healthcare environments.

Multi-surface cleaning

Can be used on both surfaces and equipment.

Convenient and effective

Disposable wipes that are ideal for general cleaning and damp-dusting.

Alcohol & disinfectant free

Safe to use on the skin with similar technology found in many shampoos and skin care products.

DETERGENT WIPES

Multi-surface wipes for general cleaning and damp dusting. They leave no smear whilst providing a superior level of cleaning power.



DETERGENT CLIP PACKS

Working with the 'Design Bugs Out' initiative, our Detergent Wipes have an integrated clip, keeping wipes neatly accessible, off the floor and nearby for whenever they are needed. The integrated clip is reinforced and angled for a firm grip and securely fits on to BP machine stands, porters chairs, bed frames, monitor stands, cabinet doors, trolleys etc.

RANGE OF USES

SURFACES

Mattresses, Bed frames, Table tops, Door handles, Desks, Chairs, Ledges, Examination couches etc.

EQUIPMENT

Telephones, Keyboards and Accessories, Wheelchairs, Trolleys, Walking frames etc.

OPTIMUM CLEANING

Clinell Detergent Wipes are unlike other detergent wipes. Our detergent formula is incredibly effective and environmentally friendly, to provide optimum, multi-surface cleaning and damp dusting in healthcare environments. Many competitor's detergent products are oil and phosphate based which leaves smears, residue and is associated with water pollution and other environmental damage.

Clinell's patent pending Detergent Wipes use compounds that lower the surface tension of a liquid and act as detergents. These are derived from renewable and sustainable plant sources such as sugar cane, rather than non-sustainable oil helping them to biodegrade rapidly.


The surfactant molecules are designed for minimum smear and excellent cleaning power.

The efficiency of these molecules also means less surfactant is necessary for a given task, contributing to minimal residue and smear.


Containing no alcohol or disinfectant, they are extremely safe to use on the skin with similar technology found in many shampoos and skin care products.



Wearing the recommended PPE, remove one wipe from the pack.

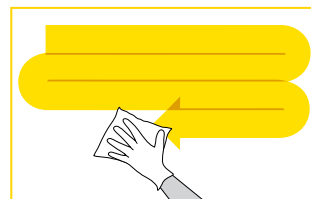

Working from clean to dirty, wipe in an 'S' shaped pattern, taking care not to go over the same area twice. Change wipe if it becomes dry or soiled and discard. Dry surface with a paper towel.

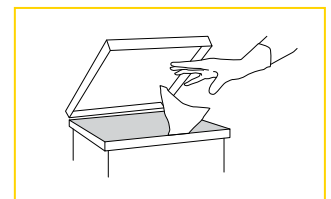

Change wipe if it becomes dry or soiled and discard. Dry the surface with a clean dry paper towel or cloth.



T: 020 7993 0030
E: info@clinell.com

| PRODUCT | UNIT OF ISSUE | ORDER CODE | NHSSC |
|-----------------------------------|---------------|------------|--------|
| Detergent Wipes Clip Pack 50 | Pack of 50 | CDCP50 | VJT254 |
| Detergent Wipes 50 (without clip) | Pack of 50 | CDW50 | VJT450 |
| Detergent Maceratable Wipes 60 | Pack of 60 | CDWM60 | VJT430 |
| Detergent Wipes 100 | Pack of 100 | CDW100 | TBC |
| Detergent Wipes Tub 100 | Tub of 100 | CDT100 | VJT225 |
| Detergent Wipes Tub Refill 100 | Refill of 100 | CDT100R | VJT226 |
| Detergent Wipes 200 | Pack of 200 | CDW200 | VJT129 |
| Detergent Wipes Bucket 225 | Bucket of 225 | CDB225 | VJT178 |
| Detergent Wipes Bucket Refill 225 | Refill of 225 | CDB225R | VJT179 |
| White Wall Mounted Dispenser | (for CDW200) | CDWDW | - |
| Wire Wall Mounted Dispenser | (for CDT100) | CDT100D | - |
| Wire Wall Mounted Dispenser | (for CDB225) | CDB225D | - |