



Fannin

Caring for life

GELTEX®inside

acute/primary care mattress

Superior performance and ultimate pressure care Acute mattress with 50,000ppm Closure chlorine resistant cover.

VERY HIGH RISK

Unique GELTEX®inside layer delivers the ultimate in patient comfort and patient care with superior pressure redistribution performance.

Applications


- Highest performance pressure care - including heel zones
- Highest patient comfort experience - breathable and supportive
- For use by any patient at VERY HIGH RISK** of developing pressure ulcers
- Demanding applications with regular cleaning requirements
- Where improved incontinence and fluid resistance is required to prevent foam core contamination
- Highest performance and weight limit requirements
- 24 hour / 7 day use
- From Domestic, Care Home and Nursing Home use to Acute and Hospital Ward environments

Features


- Layered foam construction with unique Geltex Inside top layer
- Firm edge and base support for ease of patient transfer and mattress durability
- 2 way turn to extend mattress life
- 216kg / 34 stone weight limit*
- 50,000ppm Chlorine resistant cover with welded seams, zipped closure with zip flap



with profiling interior hinges

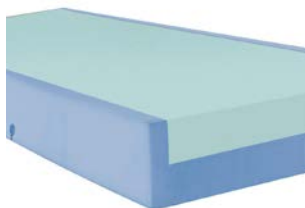

ORDERING CODE: FX-C6GT

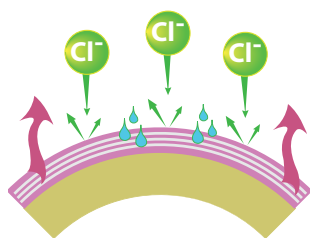

cover protection for in situ disinfection to 50,000 ppm available chlorine

GELTEX®inside
by Rectice^l®

Construction

Geltex® inside foam surface with firm edges and supportive base
Unique machine-hinged profile design, non-slip base
Complies with BS7177: 2008 medium hazard

Patient weight guide* - 216 kg / 34 stone*

Product specifics

Length - 1930mm
Width - 910mm
Depth - 150mm

Standard size - alternative sizes available. POA

Weight - 15 kg
Warranty - 4 years reducing***

Foam specification - G52/080 SB40/170

Mattress cover specification

Closure PU coated, multi-stretch, vapour permeable, waterproof knitted nylon.
High frequency welded seams and zip flap to cover zip closure.
Wipe clean cover, impregnated with anti-microbial agent in the fabric coating.
Resistant to bacterial growth, including MRSA, fungal growth and house dust mite.
White substrate fabric to help identify interior staining.

Distributed by: FANNIN LTD, Fannin House, South Business Park, Leopardstown, Dublin 18, Ireland
Tel: + 353 1 2907000 Fax: + 353 1 2907111 www.fannin.eu

* Safe working loads / weight guides are based on test methods devised at the Recticel International Development Centre, Wetteren, Belgium. Weight limit can vary on type and length of use - the weight limits are to be used as a guide only and in conjunction with clinical judgements.

** Indication of risk category - not a prescriptive criteria. *** Refer to written warranty for full details.