

Prima™ Visco mattress

Features

Closure fabric protection to 50,000 ppm chlorine

Removable cover with high frequency welded seams and zipped closure

Audit box, rotation cycle and care and cleaning instructions printed on to the cover

Very high risk specification** incorporating profiled memory foam surface

Unique machine-hinged profile design with firm edges

Rotate only design

Benefits

- Increases longevity and protects against aggressive cleaning regimes.
- Covers can be cleaned in situ for easy care or laundered.
- Best practice guidelines to aid longevity. Audit box to record inspections and rotation cycle guide information.
- Clinically tested, finest quality high resilient foams used to provide excellent pressure distribution and improve patient comfort.
- Suitable for profiling and non-profiling bed bases.
- For ease of manual handling.

Construction

Profiled memory foam surface with firm edges and supportive base
Unique machine-hinged profile design, non-slip base
Complies with BS7177: 2008 medium hazard

Patient weight guide*

- 212.9kg / 33.5 stone

Product specifics

Length - 1930mm

Width - 910mm

Depth - 150mm

Standard size – alternative sizes available

Weight - 16 kg

Warranty - 4 years reducing***

Risk category**


- Very high

Foam specification


- V50/90 / SB40/170

Mattress cover specification

Closure PU coated, multi-stretch, vapour permeable, waterproof knitted nylon. High frequency welded seams and double-sided flap to protect zip closure. Wipe clean cover, impregnated with anti-microbial agent in the fabric coating. Resistant to bacterial growth, including MRSA, fungal growth and house dust mite. White substrate fabric to help identify interior staining.

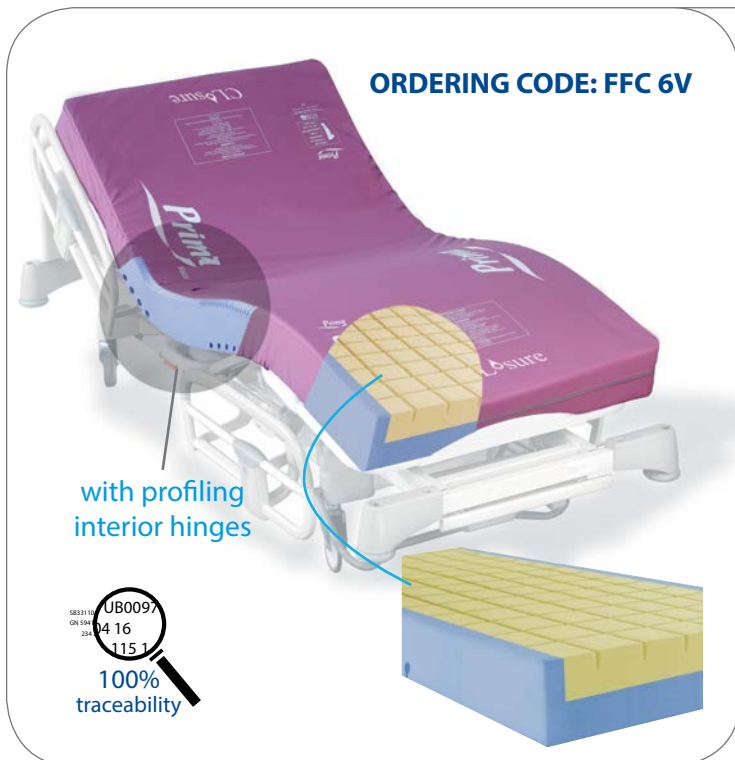

CL^osure®

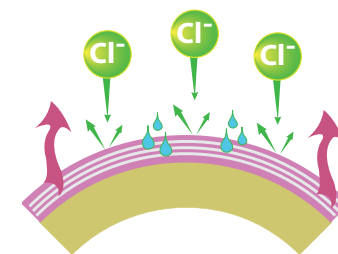

cover protection for
in situ disinfection
to 50,000 ppm
available chlorine

Detailed care & cleaning instructions are available on request.