

Absorbent Theatre Mat



Nordiska - Absorbent Theatre Mat

Suitable for use under the operating table

Dual-sided

Allows the mat to absorb, store and bind fluid through both sides

Superabsorbent Core

Absorbs and retains liquids in a safe and hygienic manner

Polypropylene non-woven fabric outer layer


Provides a fatigue resistant cover, while allowing fluid to pass through into the core

Single-use

Non-sterile

Latex-free

nordiska+

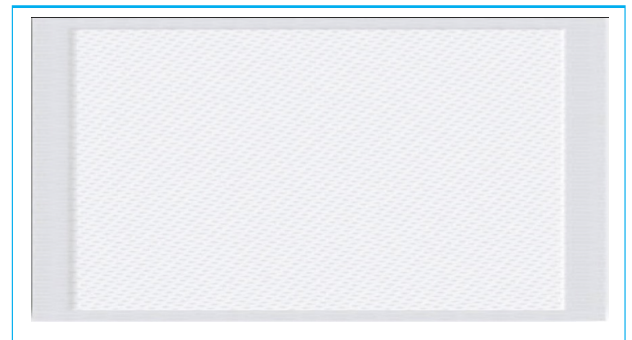

External Dimensions	75cm x 36cm
Capacity	Up to 8 litres
Pack Size	100 pieces (mats) per box
Product Code	95922
NDC SKU	235551

DCC  ITAL



Mats must be used white side up, i.e. stand on the white side, with the blue side facing down

Safehold - Anti-slip disposable absorbent floor mat - Instructions for use

The anti-slip disposable absorbent floor mat is always used with the blue (3-litre mat) side in contact with the floor, and the white side facing up... so the user is standing with their feet on the white side.

Once wet, capillary action on the blue side of the mat will hold the mat in place, preventing movement and helping to prevent slips.

This capillary action also means the mat can be thrown down on to an already wet floor - it will start to absorb and 'suck the water up'.

Once used, the mat can be safely incinerated

Dual-sided

Allows the mat to absorb, store and bind fluid through both sides

Single-use

Non-sterile

Latex-free

For more information, please contact our Fannin UK office on: 01453 874130

External Dimensions	58cm x 100cm
Capacity	3 litres per m2
Pack Size	60 pieces (mats) per box
Product Code	FHCAM2
NDC SKU	235544