



# LAPAROSCOPIC CONSUMABLES



VERESS NEEDLES



SUTURE GRASPER



DIATHERMY POUCH



DIATHERMY QUIVER

## VERESS NEEDLES

### Easi-VN Pneumoperitoneum Insufflation Needle

**Fannin UK Ltd single use pneumoperitoneum insufflation needle provides a reliable, safe, cost effective means of gas insufflation in laparoscopic surgery**

**Single Use Latex Free :** Guarantees sterility & sharpness, avoiding difficult time consuming cleaning and re-processing of reusable needles which may be blocked, blunt or malfunctioning

**Sizes:** Available in two sizes 120mm (14G) - 150mm (14G)


**Sharp Bevelled Tip:** Quick, safe and easy insertion

**Obturator :** Smooth retracting obturator minimises risk of accidental adjacent organ damage

**Injection Moulded Hub:** comfortable grip, strong and no risk of needle detachment

**Luer Lock with stop cock:** Easy connection and control of gas

**Fannin can offer both a standard and audible veress needles, the audible sound and colour indicator confirms the correct insertion of the needle.**

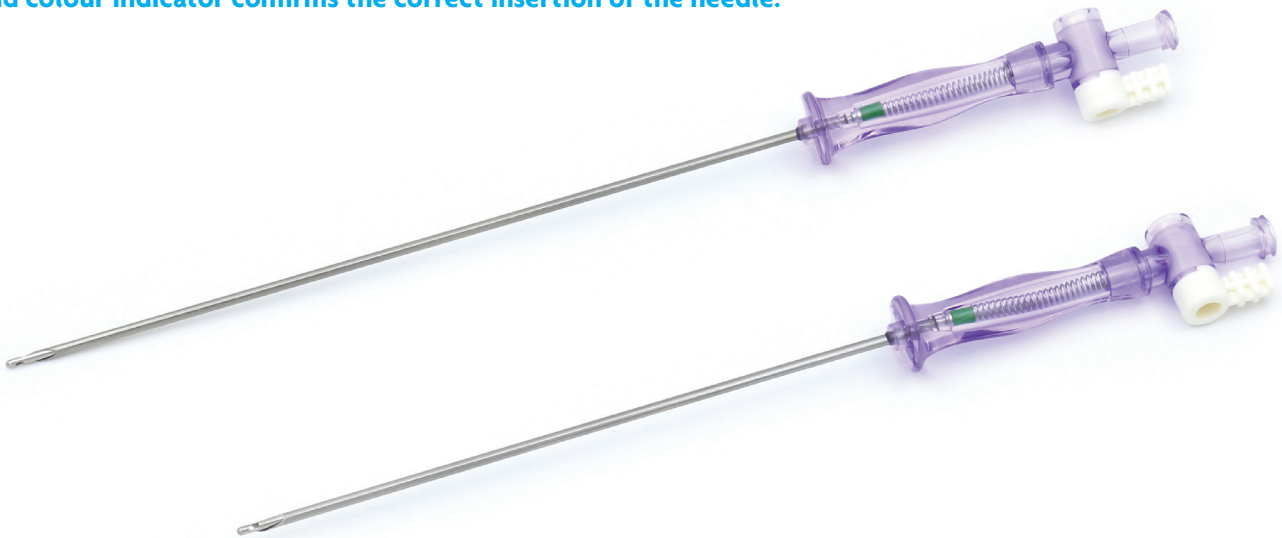

### Ordering details


Code	Description	UOI	Scotland NDC SKU Code
VN0001	Easi-VN Pneumoperitoneum Insufflation Needle 120mm (14G)	Box of 25	194360
VN0001 (NHS SC)	Easi-VN Pneumoperitoneum Insufflation Needle 120mm (14G)	Box of 50	-
VN0002	Easi-VN Pneumoperitoneum Insufflation Needle 150mm (14G)	Box of 25	194377
VN0003	Easi-VN Audible Pneumoperitoneum Insufflation Needle 120mm (14G)	Box of 50	-
VN0004	Easi-VN Audible Pneumoperitoneum Insufflation Needle 150mm (14G)	Box of 50	-

# SUTURE GRASPER

## Suture Grasper - Mediflex Surgical Products

Mediflex offers a full line of Suture Grasper devices to maximise patient safety and surgical efficiency

- 14-gauge stainless steel device attached to an ergonomic handle Innovative design securely holds suture in place
- Percutaneously pass and retrieve sutures
- One-hand function with comfortable finger platforms for tactile control Grasper arms permit easy suture access & retrieval
- Multiple clinical applications
- Accommodates 0 to 2.0 suture sizes

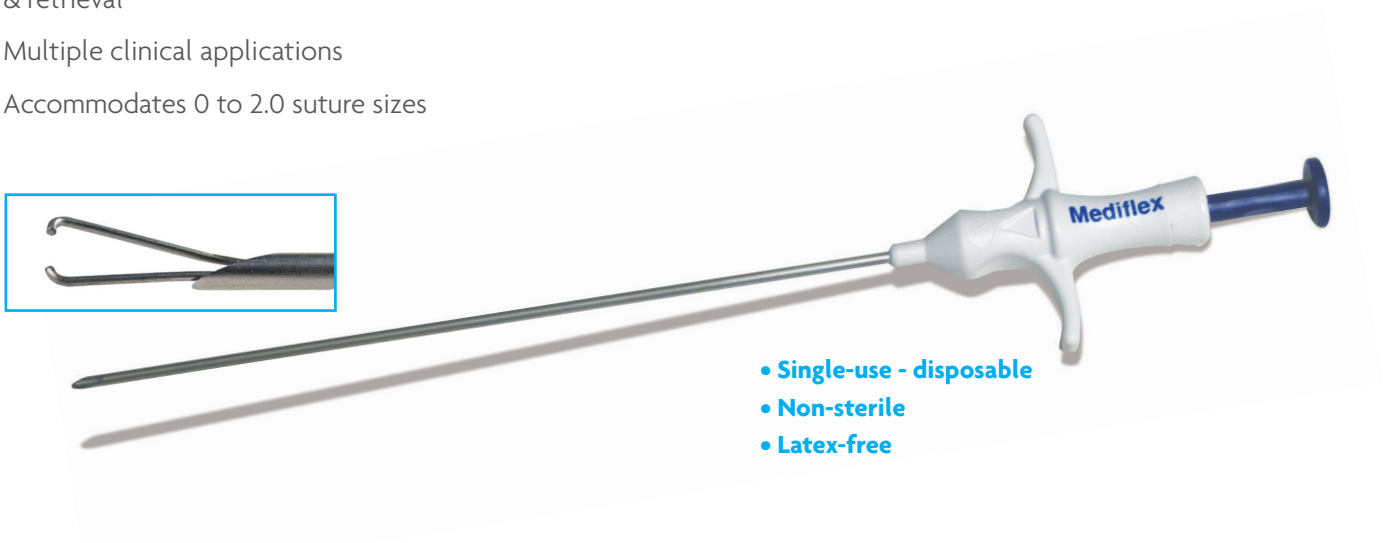

### Applications

**Port Closure** - Without extending the incision, enables the surgeon to rapidly close each laparoscopic port securely and safely while pneumoperitoneum is maintained

**Ventral Hernia Mesh Fixation** - Allows the surgeon to temporarily fixate the mesh corners or center prior to final placement; prevents mesh migration

**Organ Anchoring** - Facilitates percutaneous anchoring of organs

**Percutaneous Suturing** - Excellent choice in any procedure requiring percutaneous suturing in lieu of utilizing an additional trocar

**Vessel Ligation** - Tying a knot around a bleeding vessel too large to coagulate can provide effective haemostasis

### Ordering details

Code	Description	UOI
D1001	Mediflex Laparoscopic Suture Pass Grasper Device	Box 10

## DIATHERMY POUCH & QUIVER

### Diathermy Pouch

- **Sterile** - Single use 45cm x 33cm
- **Heavy Duty** Tough Polyethylene Instrument Pouch
- **Two** compartments
- **Strong Adhesive strip** to securely attach pouch to drape
- **Convenient:** Allows close placement of instruments in the sterile field during surgery
- **Latex Free**



### Disposable Sterile Extending Diathermy Quiver

- **Single** use
- **Double wrapped** sterile
- **Single Quiver 200mm** with optional Extension tube giving over length 400mm x 65mm
- **Lightweight**, strong with attachment loops
- **Latex Free**



## Ordering details

Code	Description	UOI	Scotland NDC SKU Code
LAPQ-0001	Laparoscopic Eextended Quiver	Box 10	219575
59149-100	Easi-Pouch Disposable sterile Instrument Pouch	Box 100	-