



LAPAROSCOPIC CONSUMABLES



VERESS NEEDLES



SUTURE GRASPER



DIATHERMY POUCH



DIATHERMY QUIVER

VERESS NEEDLES

Easi-VN Pneumoperitoneum Insufflation Needle

Fannin UK Ltd single use pneumoperitoneum insufflation needle provides a reliable, safe, cost effective means of gas insufflation in laparoscopic surgery

Single Use Latex Free : Guarantees sterility & sharpness, avoiding difficult time consuming cleaning and re-processing of reusable needles which may be blocked, blunt or malfunctioning

Sizes: Available in two sizes 120mm (14G) - 150mm (14G)


Sharp Bevelled Tip: Quick, safe and easy insertion

Obturator : Smooth retracting obturator minimises risk of accidental adjacent organ damage

Injection Moulded Hub: comfortable grip, strong and no risk of needle detachment

Luer Lock with stop cock: Easy connection and control of gas

Fannin can offer both a standard and audible veress needles, the audible sound and colour indicator confirms the correct insertion of the needle.

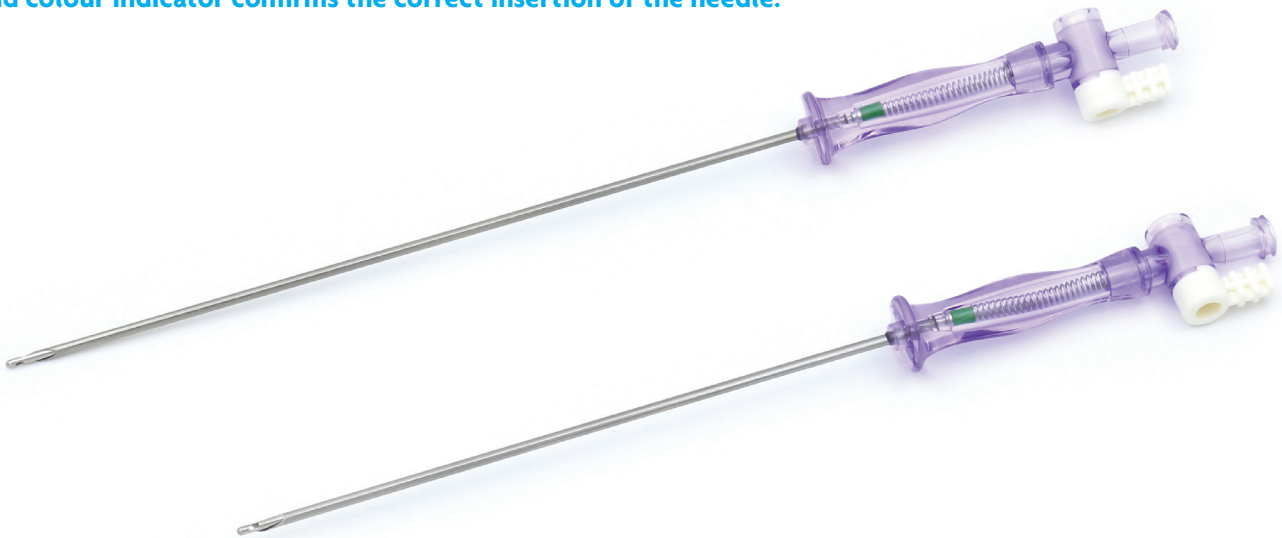

Ordering details


Code	Description	UOI	NHS Supply Chain code	Scotland NDC SKU Code
VN0001	Easi-VN Pneumoperitoneum Insufflation Needle 120mm (14G)	Box of 25	-	194360
VN0001 (NHS SC)	Easi-VN Pneumoperitoneum Insufflation Needle 120mm (14G)	Box of 50	NPC FTR 1070	-
VN0002	Easi-VN Pneumoperitoneum Insufflation Needle 150mm (14G)	Box of 25	NPC FTR 1071	194377
VN0003	Easi-VN Audible Pneumoperitoneum Insufflation Needle 120mm (14G)	Box of 50	-	-
VN0004	Easi-VN Audible Pneumoperitoneum Insufflation Needle 150mm (14G)	Box of 50	-	-

SUTURE GRASPER

Suture Grasper - Mediflex Surgical Products

Mediflex offers a full line of Suture Grasper devices to maximise patient safety and surgical efficiency

- 14-gauge stainless steel device attached to an ergonomic handle Innovative design securely holds suture in place
- Percutaneously pass and retrieve sutures
- One-hand function with comfortable finger platforms for tactile control Grasper arms permit easy suture access & retrieval
- Multiple clinical applications
- Accommodates 0 to 2.0 suture sizes

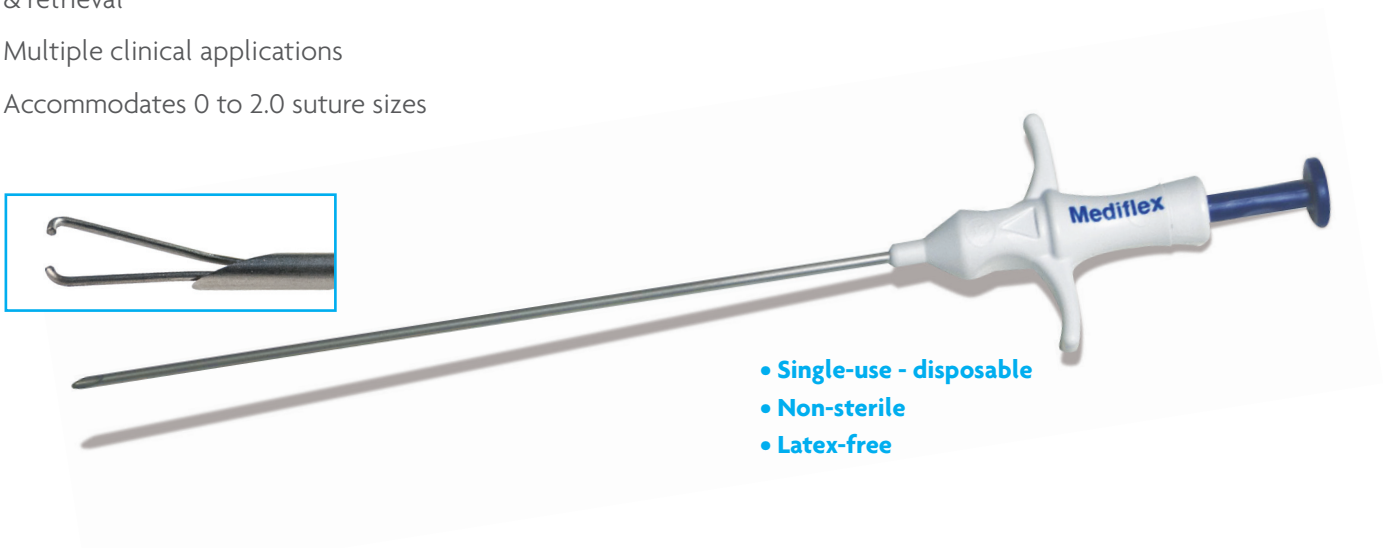

Applications

Port Closure - Without extending the incision, enables the surgeon to rapidly close each laparoscopic port securely and safely while pneumoperitoneum is maintained

Ventral Hernia Mesh Fixation - Allows the surgeon to temporarily fixate the mesh corners or center prior to final placement; prevents mesh migration

Organ Anchoring - Facilitates percutaneous anchoring of organs

Percutaneous Suturing - Excellent choice in any procedure requiring percutaneous suturing in lieu of utilizing an additional trocar

Vessel Ligation - Tying a knot around a bleeding vessel too large to coagulate can provide effective haemostasis

Ordering details

Code	Description	UOI
D1001	Mediflex Laparoscopic Suture Pass Grasper Device	Box 10

DIATHERMY POUCH & QUIVER

Diathermy Pouch

- **Sterile** - Single use 45cm x 33cm
- **Heavy Duty** Tough Polyethylene Instrument Pouch
- **Two** compartments
- **Strong Adhesive strip** to securely attach pouch to drape
- **Convenient:** Allows close placement of instruments in the sterile field during surgery
- **Latex Free**



Disposable Sterile Extending Diathermy Quiver

- **Single** use
- **Double wrapped** sterile
- **Single Quiver 200mm** with optional Extension tube giving over length 400mm x 65mm
- **Lightweight**, strong with attachment loops
- **Latex Free**



Ordering details

Code	Description	UOI	NHS Supply Chain code	Scotland NDC SKU Code
LAPQ-0001	Laparoscopic Eextended Quiver	Box 10	-	219575
59149-100	Easi-Pouch Disposable sterile Instrument Pouch	Box 100	NPC VJD3560	-