

Smoke Evacuation System & Accessories



Features:

- Three procedure operating modes (open tubing, laparoscopic, Electro-surgical)
- Fifty hours of filter life
- Touch screen based
- Multi-port design
- Quietest Performance
- Includes Foot Switch
- EvaQMax Auto (Automatic Activation Device) is an optional feature
- Fluid Trap is also an optional feature

EvaQMax

Smoke Free Evacuation


Procedure Operating Modes

The desired procedure can be set simply by touching the corresponding image. Automatically customizing suction default level based on desired surgical procedure.



Longer lasting filter life (up to 50 Hours)

Can save on running costs.

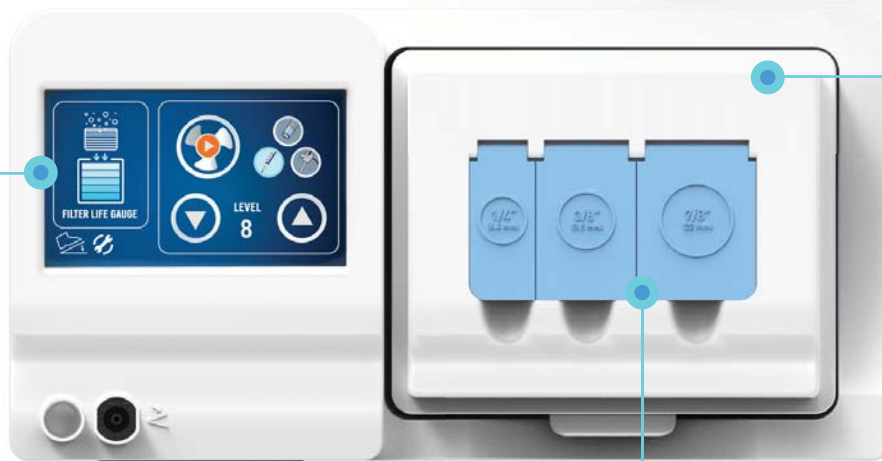


Touch screen

The touch screen panel can be operated while wearing gloves.



Multi Port design

Three different sized connection ports are designed to create an air-tight bond with the desired port. The port cover cap eliminates any possible contamination with the glove when connecting the hose to the unit.

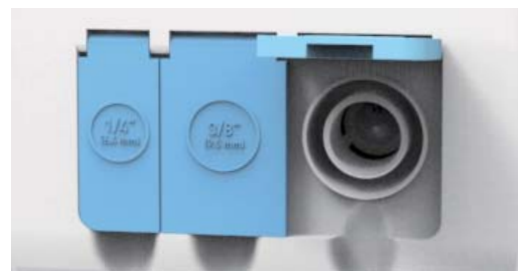

Technical Specification

Item	Specification
Noise Levels (Nominal Values)	maximum power < 55 dBA
Dimensions (H x W x D)	6 in x 11 in x 15.5 in (15.2 cm x 27.9 cm x 39.4 cm)
Flow Rates	36 cfm (708 lpm) using 7/8 in (22 mm) tubing
Weight	With filter installed: 12.0 lbs (5.0 kg) Without filter installed: 9.5 lbs (4.0 kg)
Certifications	CE, FDA
Warranty	Three (3) Year Full Warranty

Accessories



Footswitch (Included)



Automatic Activation Device (Optional)



Fluid Trap (Optional)




Model	Product Description	Packaging Info.
EQM-A	Automatic Activation Device (AAD)	1 pcs/Box
EQM- FT	Fluid Trap	20 pcs/Box
EQM-CAF	Filter for EQM-U	4 pcs/Box

Telescopic Smoke Evacuation Pencil

ONE TOUCH BUTTON EXTENDABLE TELESCOPIC TYPE

The easiest way to extend the length

1. Ability to keep focus on the surgery area even when extending the length.
2. Easy extension and retraction of suction tube with blade attachment.
3. Thin, light weight design provides comfortable usage for surgeons.
4. Superior suction Capability.

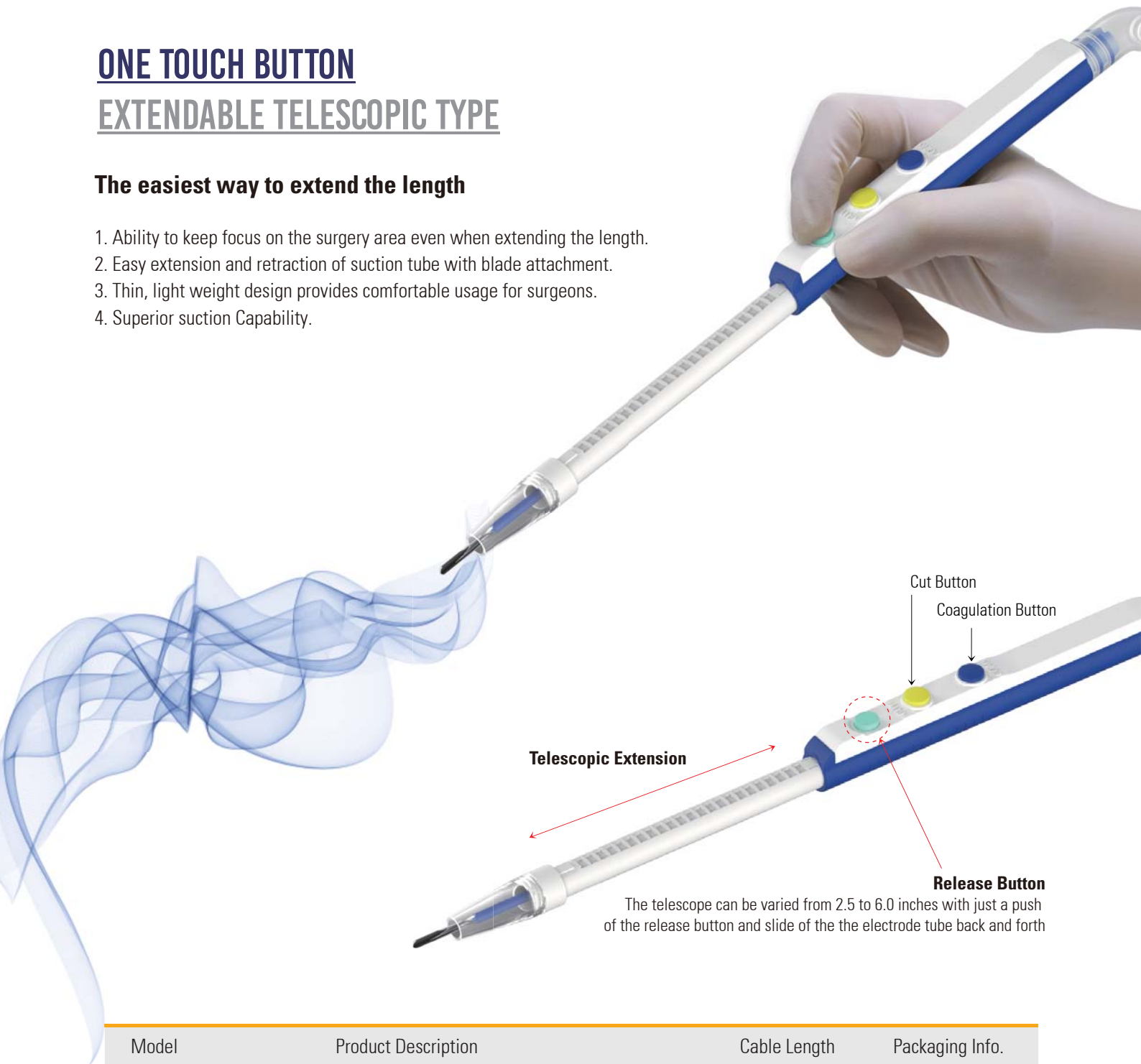

Model	Product Description	Cable Length	Packaging Info.
EvBNH	Telescopic, 2.5" Non-Stick Ceramic Coated Blade, Holster	10' Cable	1pc / pack 10 pcs/ Box 4 Boxes / Case

Dual Extension Smoke Evacuation Pencil



ULTIMATE SOLUTION for SMOKE FREE in OR DUAL EXTENSION SMOKE EVACUATION PENCIL

1. EASY TO EXTEND THE TELESCOPIC TUBE WITH THE RELEASE BUTTON

The telescopic extendable tube is convenient and easy to adjust to desired length.




2. EXTENSION OF THE SUCTION TIP ALLOWS MAXIMUM SMOKE EVACUATION

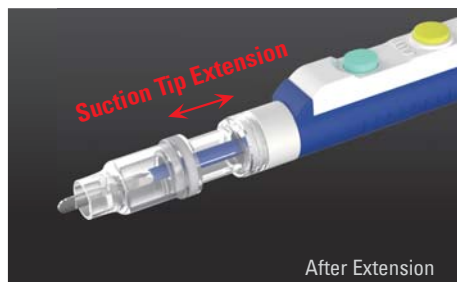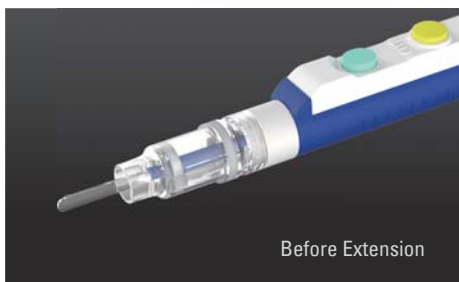



Hold the suction tip and slide it as close as possible to the blade in order to maximize the suction capacity.

Model	Product Description	Cable Length	Packaging Info.
EvDE	Dual Extension, Telescopic, 2.5" Non-Stick Ceramic Coated Blade, Holster	10' Cable	1pc / pack 10 pcs/ Box 4 Boxes / Case

Surgical Smoke Evacuation Adapter

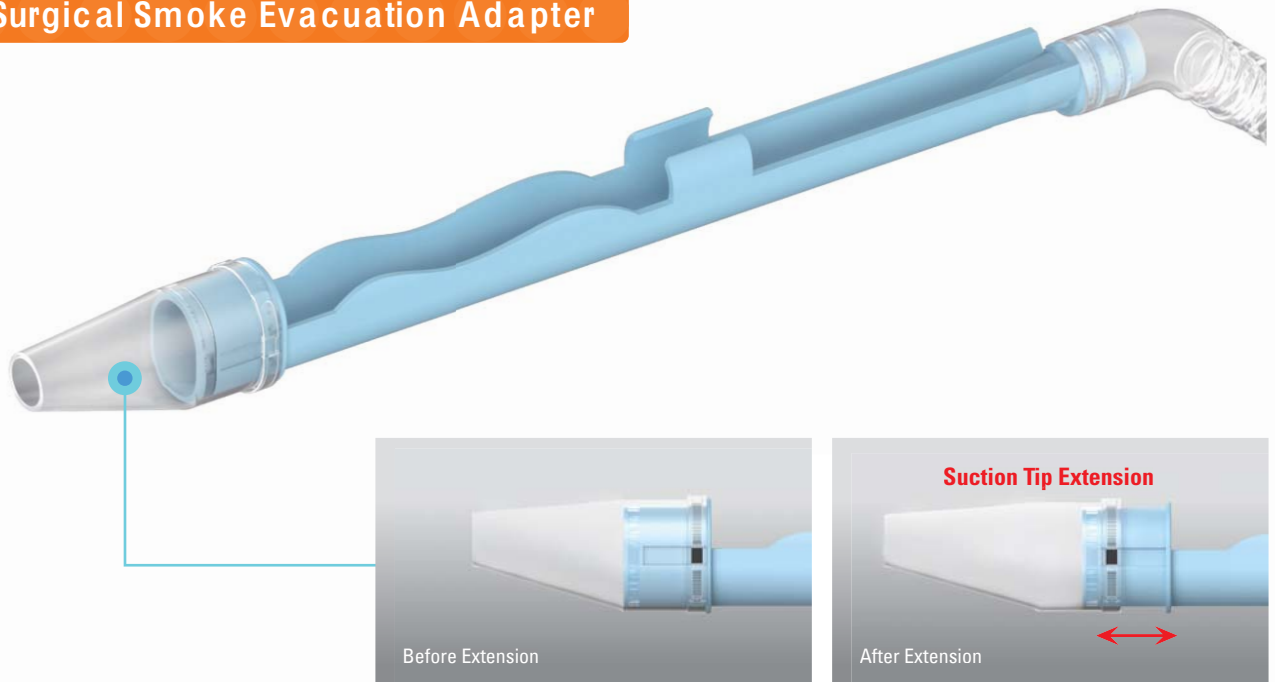

EXTENSION OF THE SUCTION TIP ALLOWS MAXIMUM SMOKE EVACUATION

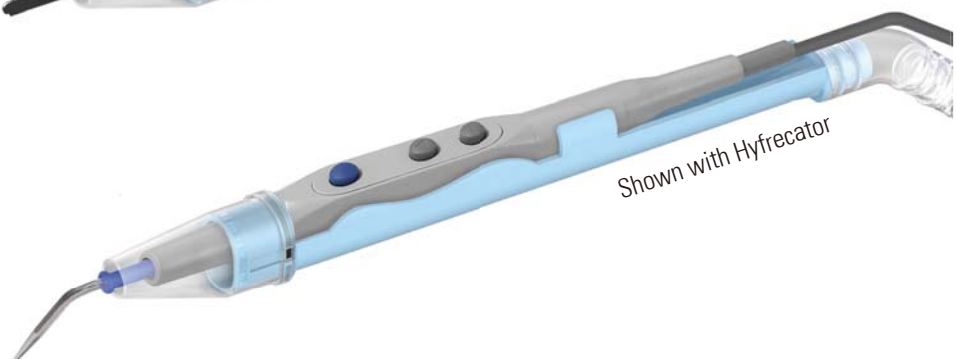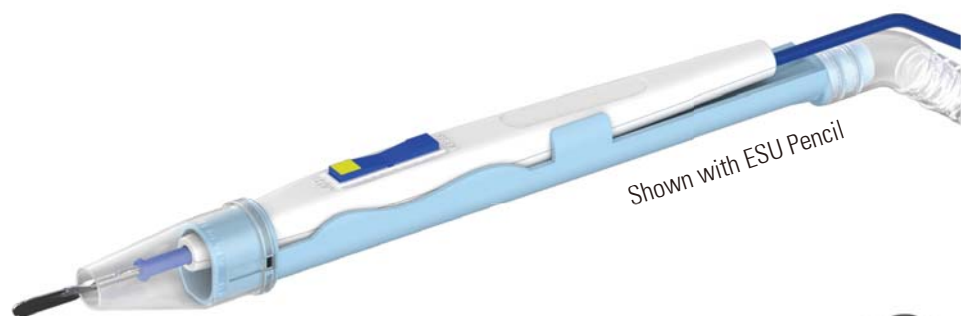

Model	Product Description	Tubing Length
EvFIT	Ideal for surgical smoke evacuator applications, Extension	3/8in X 10' Ft Tubing

Standard Smoke Evacuation Pencil - PS07

Compatible with any electrode on the market, including hex-shaped electrodes. Compact design and low-profile front end enhances visibility of the surgical field. Transparent nozzle minimizes any obstruction of the blade's tip or target tissue. Nozzle X-Ray identifiable

Cord set housed within the tubing to eliminate clutter at the surgical site

Large buttons with excellent Hand-feel

Standard
70mm / 2.75in

Angled connection and Ultra-flexible corrugated tubing optimizes flexibility and range of motion.

Silicone subtle texture enhances grip.

Slim profile for ergonomic comfort.
Significantly smaller than competitive smoke pencils.
Aerodynamic design with rapid smoke flow and less noise.

Model	Product Description	Box Quantity
PS07-K70GY	70mm Non-stick Blade electrode,S nozzle 3.0m(10 ft) Tubing, 22mm connector, 3.0m(10 ft) Cord	20PCS/Box

Telescopic Smoke Evacuation Pencil - PS11

Position the smoke nozzle anywhere without changing of the electrode

No additional lengths of electrodes or different size smoke evacuation nozzles are needed.

Square locking mechanism prevents electrode from moving during operation

Electrode is significantly more stable at maximal length position than competitive products

Compact design and low-profile front end enhances visibility of the surgical field.

Transparent nozzle minimizes any obstruction of the blade's tip or target tissue.

Angled connection and Ultra-flexible corrugated tubing optimizes flexibility and range of motion.

Large buttons with excellent Hand-feel

Cord set housed within the tubing to eliminate clutter at the surgical site

Subtle texture enhances grip.

100mm

22mm smoke evacuator connector

3-pin plug


Model	Product Description	Box Quantity
PS11-K50GY	50mm Non-stick blade electrode, Telescopic Function with locking mechanism 3.0m(10 ft) Tubing, 22mm connector, 3.0m(10 ft) Cord	20PCS/Box

Smoke Evacuation Accessories

Smoke Evacuation Liquid Trap

Catalogue # P-FT1

- 1) 22mm male inlet for Smoke Evacuation Pencil
- 2) Liquid Trap
- 3) Sponge liquid collector (50ml)
- 4) 22mm female outlet for Smoke Evacuator

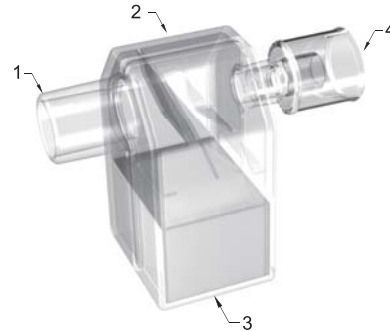

Surgical Smoke Evacuation Tubing

Catalogue # ESE-TS4

Laparoscopic, 3 meters,
 With regulator
 Sterile
 20PCS/Box



Male Luer connection



Suction regulator from 0%-100%
 connection



7/8 in (22 mm) female
 connection

Passive Smoke Evacuation Filter

Catalogue # PSF0004

NHS SC PVC Code FGB087 25PCS/Box



- Cost effective solution
- Sterile - Single use - Easy peel pouch
- Easy to use: Attached to the trocar by universal leak free luer connector
- Activated Charcoal / Bacterial/Viral Filter - 99.999% efficiency
- Thumb Wheel control
- Latex free
- Recommended for use during minimally invasive surgery involving insufflation, electrocautery, ultrasonic scalpel & laser

Fannin Ltd.

Fannin House, South County Business Park
 Leopardstown, D18 Y0C9
 T: +353 (0)1 2907000
 E: ordersl@fannin.eu
 W: www.fannin.eu

Fannin UK Ltd.

Fernhouse - Suite 3 Daniels Industrial Estate
 Bath Road, Gloucestershire
 GL5 3TJ, United Kingdom
 T: +44 (0) 1453 874 130
 E: salesll@fannin.eu