

# Master E-Sac

The fully automated  
Tissue Retrieval Sac

## DEPLOY

- using Espiner's unique introducer system

## POSITION

- using ergonomic handle

## RETRIEVE

- via laparoscopic incision reducing hospital stay, cost and complication



Safe and secure  
closure with strong  
monofilament  
drawstring



One-piece tear  
resistant design with  
unique PUR coating  
for maximum strength  
and durability




Unique  
introducer  
system



# Master E-Sac Features:

- Comprehensive solution covering all laparoscopic procedures including large spleen, ovary and kidney.
- Minimally invasive; reducing length of stay, financial strain and improved patient experience.
- Our unique coated ripstop nylon is resistant to tearing if punctured.
- Includes X-ray opaque marker.
- Polyurethane coating impermeable to biological fluids.
- Safe closure with strong monofilament drawstring.
- Introducer system and arms permit external manipulation of the Sac for tissue capture.
- Ergonomic and non-slip introducer system can be withdrawn from patient, leaving tissue securely contained until retrieval.
- Up to a five-year shelf life.

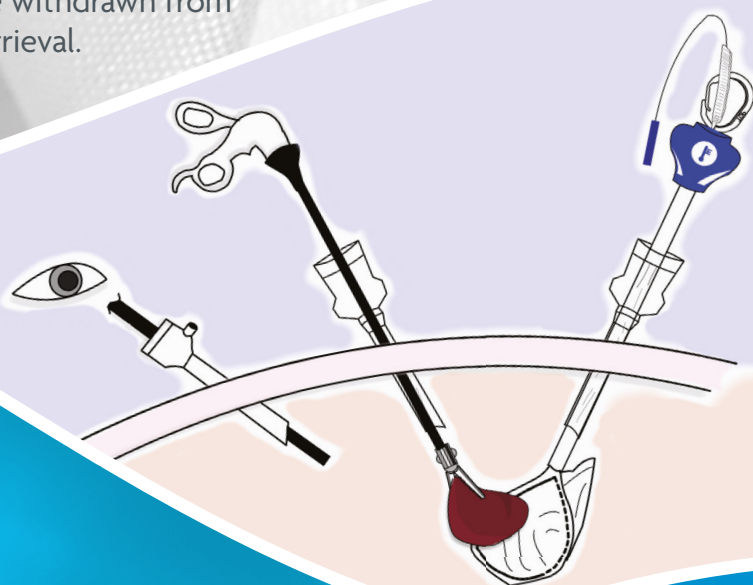

Product Code	Trocar (mm)	Diameter (mm)	Length (mm)	Volume (ml)	Number Per Box
EMP80-WG-6	10	64	105	150	6
EMP100-WG-6	10	64	150	200	6
EMP120-WG-6	10	85	210	550	6
EMP160-WG-3	12	110	260	1500	3

To talk to Espiner about our unique Espiner Master E-Sac range, please call us on **+44 (0)1275 878 801** or email **[espiner@fannin.eu](mailto:espiner@fannin.eu)**  
**[www.espinarmedical.com](http://www.espinarmedical.com)**