

EvaQMax

Smoke Free Evacuation




Procedure operating modes

The desired procedure can be selected simply by touching the corresponding image, automatically customising suction default level based on desired surgical procedure.

Longer lasting filter life

Up to 50 hours – can save on running costs.

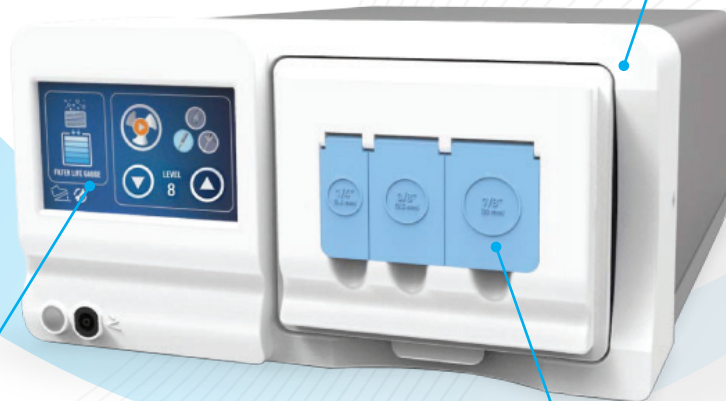


Touchscreen

The touchscreen panel can be operated while wearing gloves.



Multi-port design

Three different sized connection ports are designed to create an airtight bond with the desired port. The port cap eliminates any possible contamination with the glove when connecting to the unit.

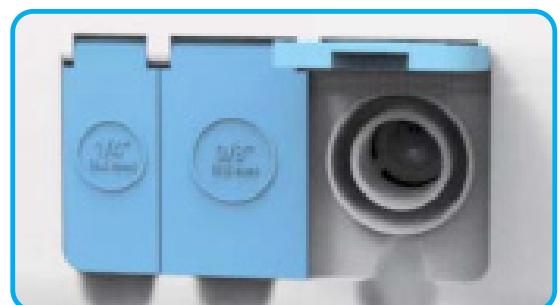

Technical Specifications

Item	Specification
Noise Levels (Nominal Values)	Maximum power <55dBA
Dimensions (H x W x D)	15.2 cm x 27.9 cm x 39.4 cm
Flow Rates	25 cfm (708 lpm) using 22mm tubing
Weight	With filter installed: 5kg Without filter installed: 4kg
Certifications	CE, FDA
Warranty	Three year warranty



Footswitch (included)



Automatic Activation
Device (optional)



Fluid trap (optional)




MPC Code	Product Description	Quantity	NPC Code
FC-U	Care E-Vac Smoke Evacuation Machine	1 per box	FBU89100
FC-A	Care E-Vac Automatic Activation Device	1 per box	FBU89102
FC-FT	Fluid Trap	20 per box	FCN89125
FC-CAF	Care E-Vac Smoke Evacuation Machine Filter	4 per box	FBU89101

Telescopic Smoke Evacuation Pencil

One touch button extendable telescopic type

The easiest way to extend the length

1. Ability to keep focus on the surgery area even when extending the length
2. Easy extension and retraction of suction tube with blade attachment
3. Thin, lightweight design provides comfortable usage for surgeons
4. Superior suction capability

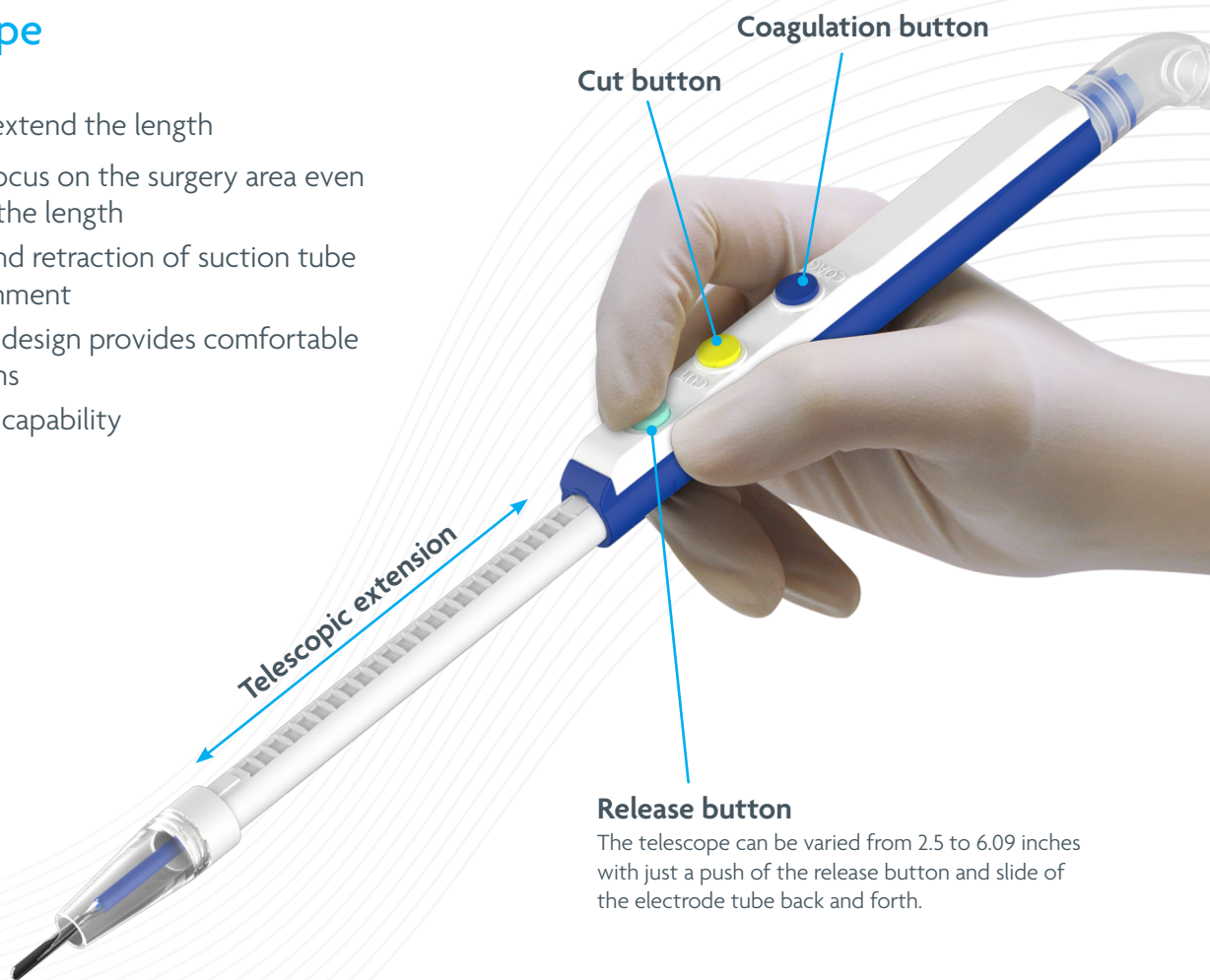

Release button

The telescope can be varied from 2.5 to 6.09 inches with just a push of the release button and slide of the electrode tube back and forth.

MPC Code	Product Description	Quantity	NPC Code
FC-BNH	Care E-Vac Telescopic Smoke Evacuation Pencil, 50 mm Non-Stick Ceramic Blade, Telescopic Function with Button Lock Mechanism, 3m (10ft) Tubing, 22mm Connector, 3m (10ft) Cord	1pc/pack 20pcs/box 4 boxes/case	FCN89127

Dual Extension Smoke Evacuation Pencil

Ultimate solution for smoke-free during surgery

Easy to extend the telescopic tube with the release button

The telescope can be varied from 2.5 to 6.09 inches with just a push of the release button and slide of the electrode tube back and forth.




Extension of the suction tip allows maximum smoke evacuation & maximum visibility

Hold the suction tip and slide it as close as possible to the blade in order to maximise the suction capacity.

MPC Code	Product Description	Cable Length	Quantity	NPC Code
FC-DE	Care E-Vac Dual Extension Smoke Evacuation Pencil 50mm, Non-Stick Blade Electrode, Telescopic Function	10" cable	1pc/pack 20pcs/box 4 boxes/case	Available direct

Surgical Smoke Evacuation Adapter


Practical solution for smoke-free surgery using your existing electro-surgical pencils

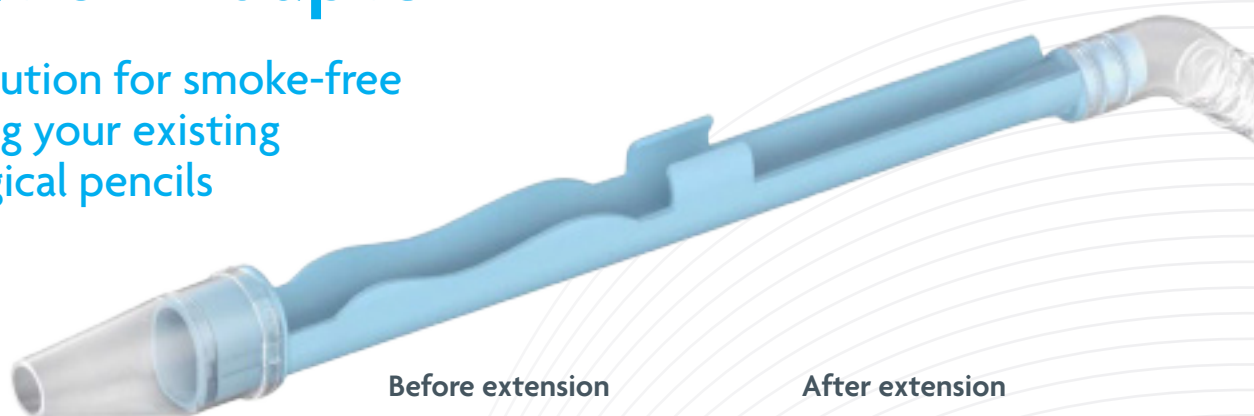

Before extension

After extension


Extension of the suction tip allows maximum smoke evacuation

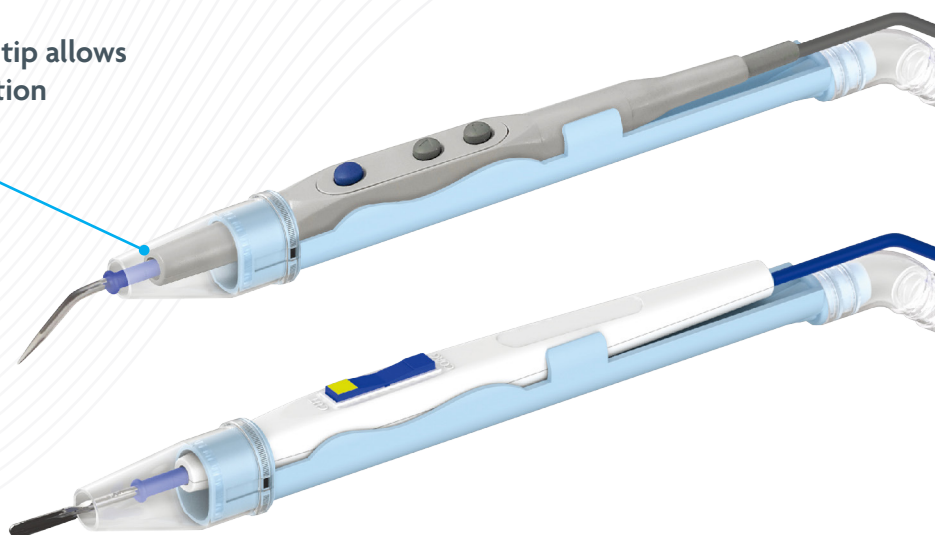

MPC Code	Product Description	Tubing Length	NPC Code
FC-FIT	Cradle Adaptor – Allows use of Standard Electro-Surgical Pencils	3/8 in x 10 ft	FBU89093

Standard Smoke Evacuation Pencil

- Compatible with any electrode on the market, including hex-shaped electrodes
- Compact design and low-profile frontend enhances visibility of the surgical field
- Transparent nozzle maximises visibility of the blade's tip or target tissues
- Nozzle x-ray identifiable tip
- Slim profile for ergonomic comfort
- Aerodynamic design with rapid smoke flow and less noise

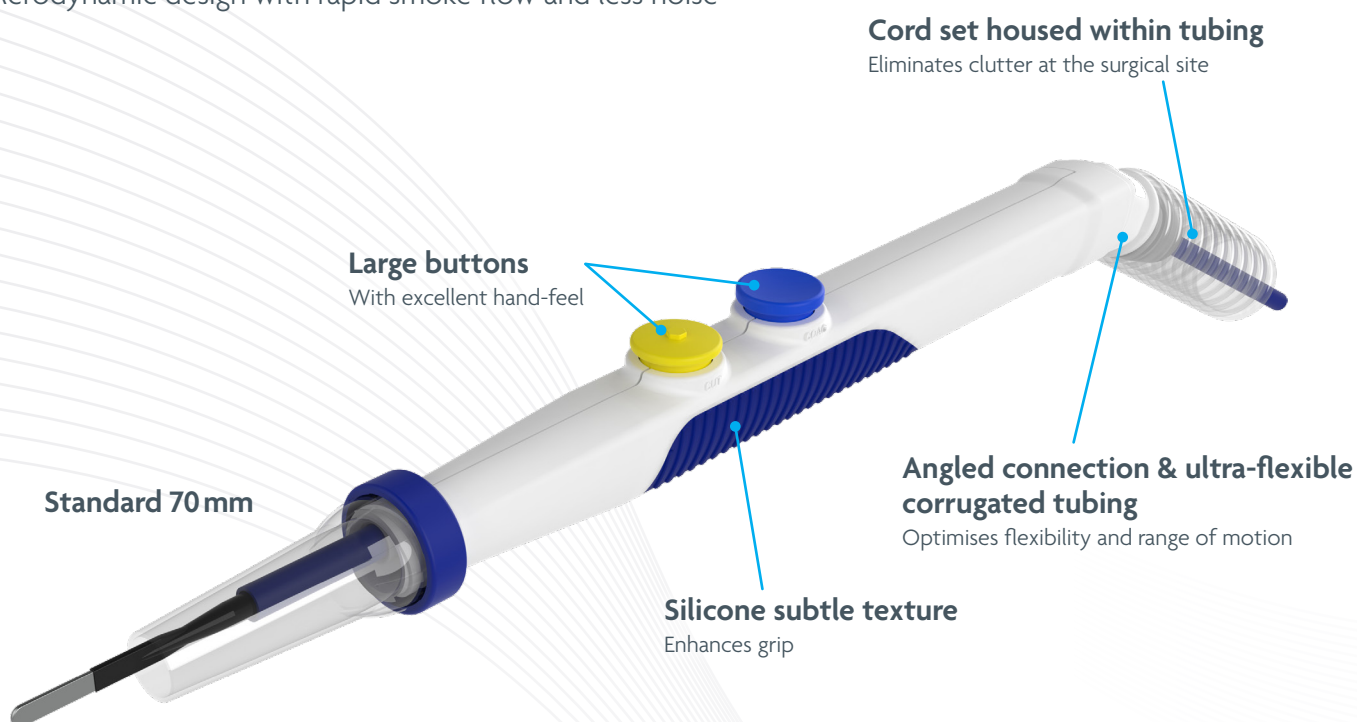

MPC Code	Product Description	Box Quantity	NPC Code
FC-PS07-K70GY	Care E-Vac Standard Smoke Evacuation Pencil 70 mm Non-Stick Blade Electrode, S Nozzle 3 m (10 ft) Tubing, 22 mm Connector, 3 m (10 ft) Cord	20 pcs/box	FBU89094

Telescopic Smoke Evacuation Pencil

- Position the smoke nozzle anywhere without changing the electrode
- No additional lengths of electrodes or different size smoke evacuation nozzles are needed
- Square locking mechanism prevents electrode from moving during operation
- Electrode is significantly more stable at maximum length position than competitor products

Angled connection & ultra-flexible corrugated tubing
Optimises flexibility and range of motion

Large buttons
With excellent hand-feel

Compact design & low-profile front-end
Enhances visibility of the surgical field

Cord set housed within tubing
Eliminates clutter at the surgical site

Transparent nozzle
Maximises visibility of the blade's tip or target tissue

Subtle texture
Enhances grip

22mm smoke evacuator connector


3-pin plug

MPC Code	Product Description	Box Quantity	NPC Code
FC-PS011-K50GY	Care E-Vac Telescopic Smoke Evacuation Pencil, 50mm Non-Stick Blade Electrode, Telescopic Function with Twist Lock Mechanism, 3m (10ft) Tubing, 3m (10ft) Cord	20pcs/box	FBU89096

Smoke Evacuation Accessories




Smoke Evacuation Liquid Trap

1. 22mm male inlet for smoke evacuation pencil
2. Liquid trap
3. Sponge liquid collector (50ml)
4. 22mm female outlet for smoke evacuator

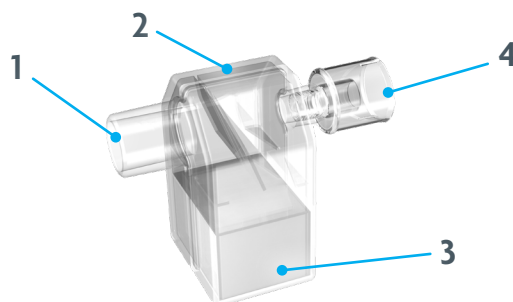

Surgical Smoke Evacuation Tubing

- Laparoscopic, 3m, with regulator
- Sterile
- 20 pieces per box

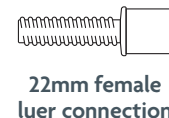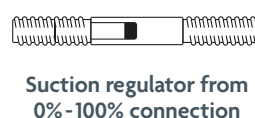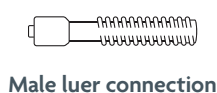

Passive Smoke Evacuation Filter

- Cost-effective solution
- Sterile, single-use, in easy-peel pouch, latex-free
- Easy-to-use – attached to the trocar by universal leak-free luer connector
- Activated charcoal/bacterial/viral filter – 99.999% efficiency
- Thumb wheel control
- Recommended for use during minimally invasive surgery involving insufflation, electrocautery, ultrasonic scalpel or laser



MPC Code	Product Description	NPC Code
FC-P-FT1	Care E-Vac Smoke Evacuation Fluid Trap 22mm 50ml	FBU89095
FC-P-F2	Care E-Vac Smoke Evacuation HEPA Filter & Fluid Trap 22mm 50ml	
FC-ESE-TS1	Care E-Vac Gynaecological Smoke Evacuation Tubing with HEPA filter 22mm connector 3m - sterile	FBU89091
FC-ESE-TS4	Care E-Vac Laparoscopic Surgical Smoke Evacuation Tubing 22mm connector 3m, with Roller Clamp Control - sterile	FBU89092
FC-ESE-TS6	Care E-Vac Open Surgical Smoke Evacuation Tubing 22mm connector 3m, with Mesh Guard - sterile	FCN 89129
PSF0004	Passive Smoke Filter	FGB087

Scan to view our
online catalogue:



Fannin Ltd
Fernhouse - Suite 3
Daniels Industrial Estate
Bath Road, Gloucestershire
GL5 3TJ, United Kingdom
T: +44 (0) 1453 874 130
E: serviceuk@fannin.eu
fannin.eu