



FLOWArt[®] One Directional Needle-Free Valve

Neutral displacement design with integrated
check valve & silicone split septum

- One-directional flow
- Clear transparent housing to allow fluid pathway clarity and enhanced flushing
- Reduces heparin flush requirement
- Avoids blood clotting in IV catheter
- Minimises blood aspiration
- Prevents backflow during IV therapy caused by patient movement and pressure differences
- Conforms to EN ISO 80369-7
- DEPP-free, PVC-free, and latex-free


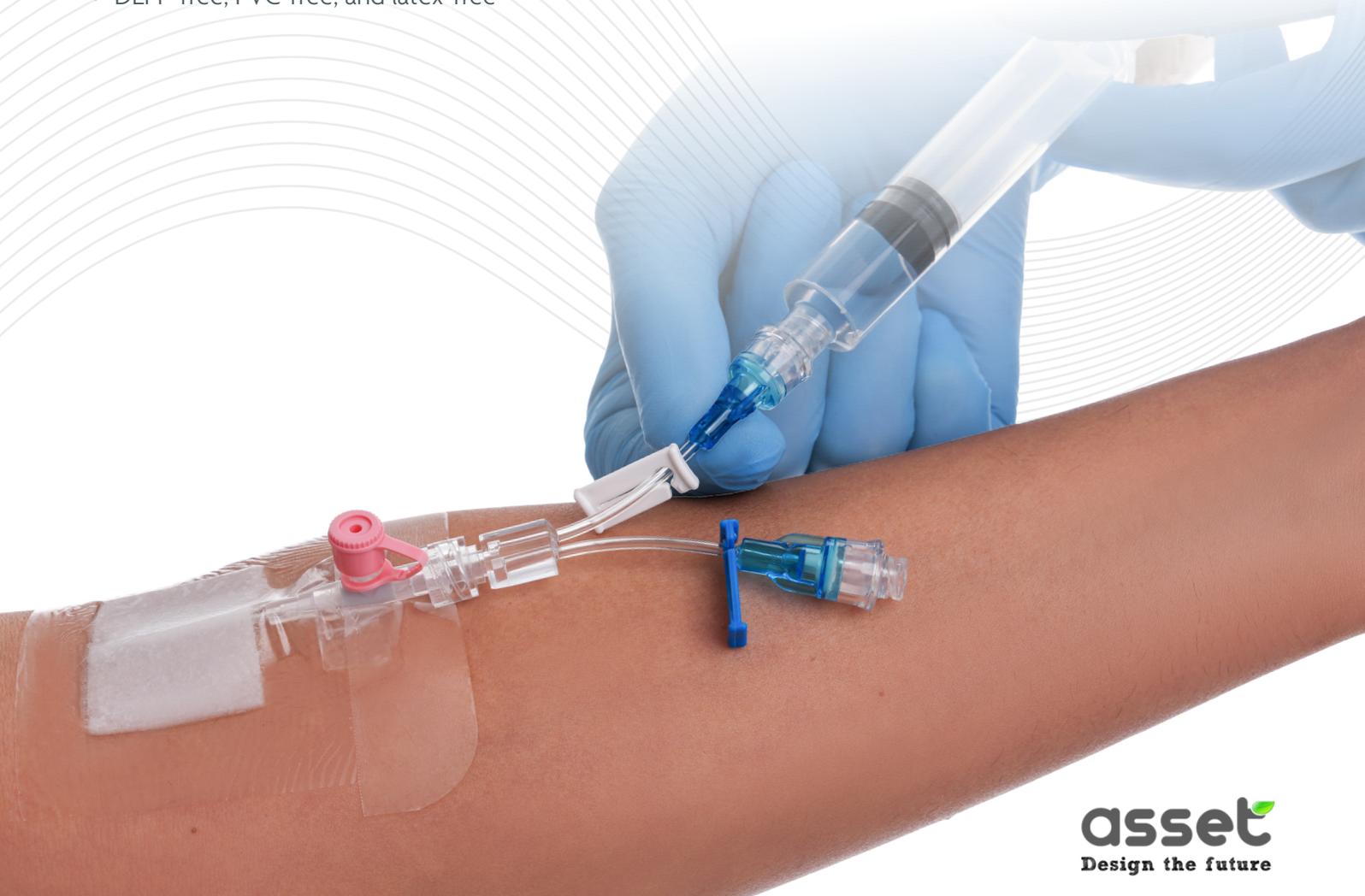


FlowArt® One Directional Needle-Free Valve has an integrated flat silicone seal that prevents blood back flow while achieving one directional flow

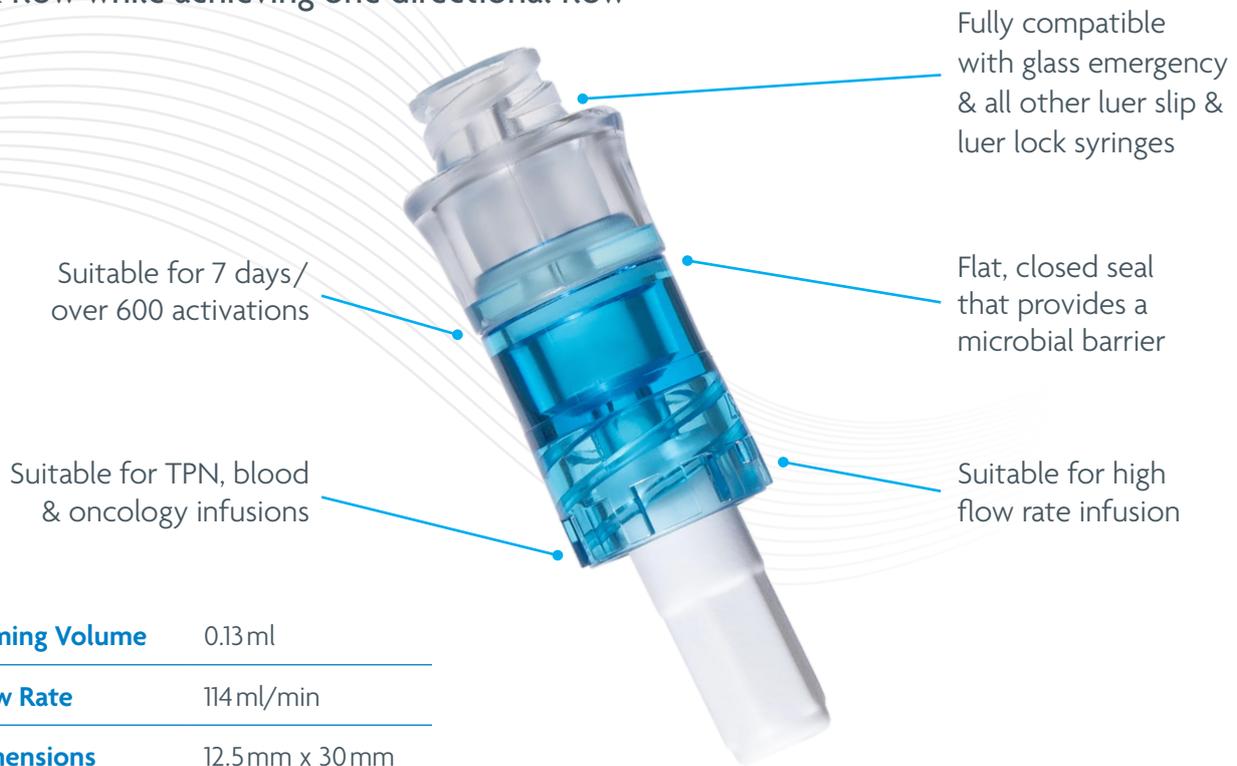

Priming Volume	0.13ml
Flow Rate	114 ml/min
Dimensions	12.5mm x 30mm
Displacement	0.01ml (neutral)

Ordering Details

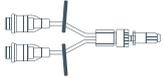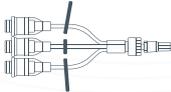Product Code	Design	Product Description	Length (cm)	NHSSC Code
A1010TY		One Directional Needle-Free Valve	3cm	—
AU1010TY		One Directional Needle-Free Valve with Extension Line	8cm	FNY85059
A2010TY		Double Lumen One Directional Needle-Free Valve	8.5cm	FNY85081
A3010TY		Triple Lumen One Directional Needle-Free Valve	8.5cm	FNY85063

