
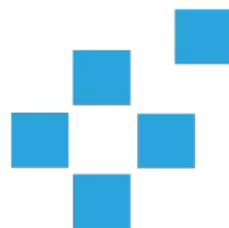


TECNO DRAPE®

ANTIMICROBIAL SURGICAL INCISE DRAPE


TecnoDrape is effective and essential protection against the contamination of the operating field and surgical incision site. TecnoDrape is designed to prevent the transmission of pathogens from the non-sterile skin towards the sterile surgical field in order to reduce the risk of the patient surgical site infections.

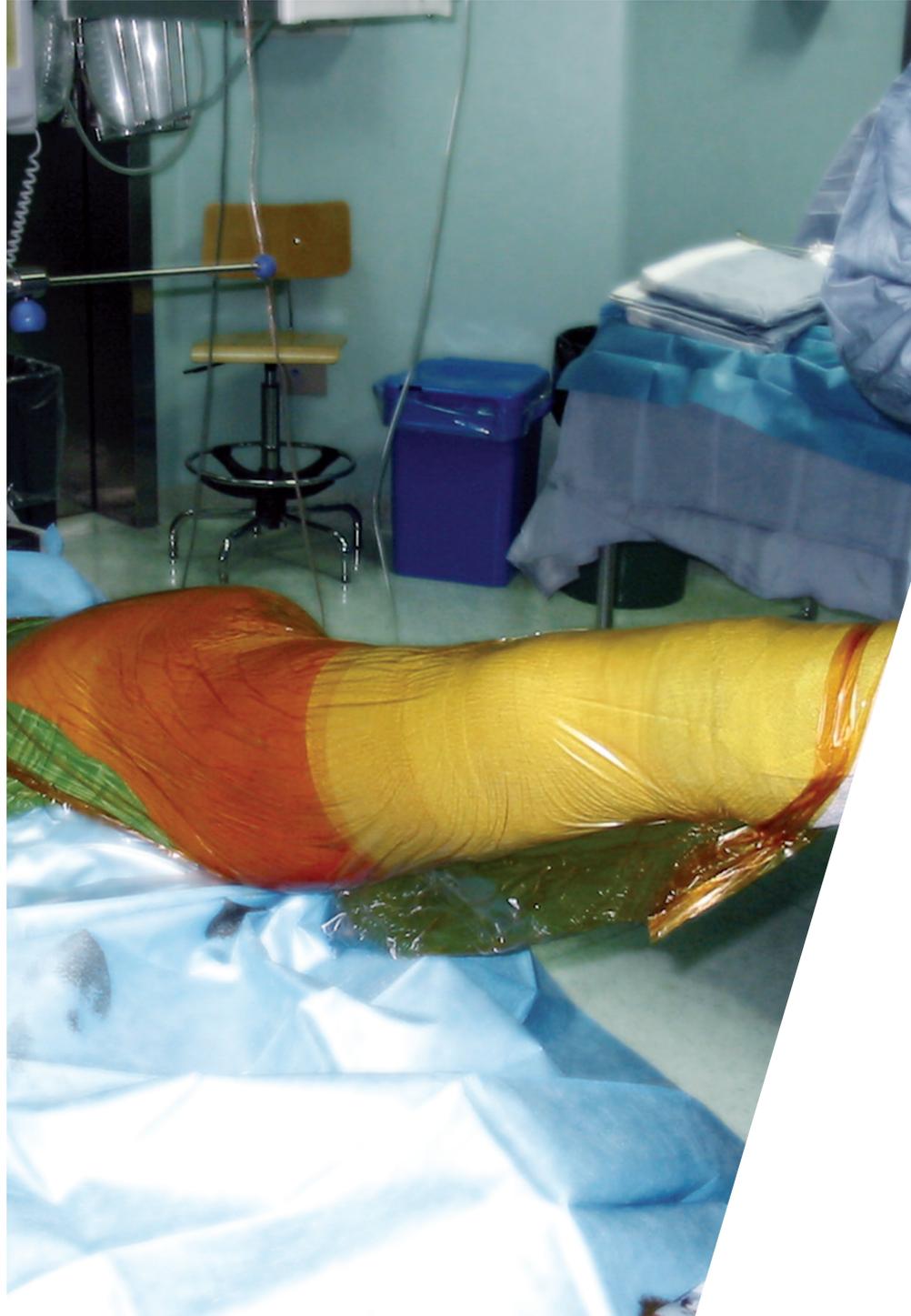

Fannin

TECNODRAPE is a surgical incise drape made with **POVIDONE IODINE** and it has a grip tape made of a slow release iodine complex. **TECNODRAPE** is effective protection against the contamination of the operative field and the surgical incision site.

FEATURES	BENEFITS
Iodine Solution	Broad spectrum bactericidal action. Fights deep bacterial flora into the top cell layers of the skin.
Gradual release of Iodine particles	Sustained duration of action with the gradual release of iodine molecules for up to six hours from the time of first application.
Acrylate adhesive tape	Strong adhesion. Perfect adhesion to skin even in very long surgeries. Perfectly isolates the cut borders preventing microorganism's migration. Lock and control of other drapes. No skin damage after removal. Hypoallergenic. Excellent tolerance for the skin.
Polyurethane support	High conformability. Emo-hydro-alcohol repellent. Impenetrable barrier against bacteria. Antiglare. RX Transparent. High extensibility (600%). No memory. Antistatic non-conductor.

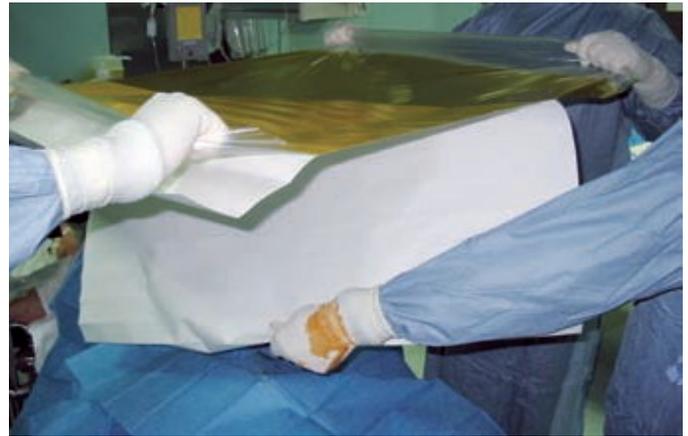

PACKAGING

- Single pack with double sterile barrier system
- Easy-opening to avoid cross-contamination
- 2 removable labels for easy traceability
- Internal dispenser for clean storage
- Dimensions and weight are certificated according with Dlgs 81/2008 ex626/94

TECNODRAPE MEDICAL CLASSIFICATION

TECNODRAPE is classified as a Medical Device CE Class III according to directive 93/42/ CEE, notifying body CE0373.

Class III All the devices include a substance that could be considered a medical drug if used separately and by itself, according to article 1 of 65/65/CEE, and that could cause incidental effects to the human body that fall within the III class.

NHS SC Code	Product Code	Total Drape Size	Incise Area	Units per Case
VEB961	IH6635TD	15 x 20cm	10 x 20cm	50
VEB959	IH6640TD	60 x 35cm	35 x 35cm	50
VEB957	IH6648TD	90 x 60cm	60 x 60cm	50
VEB960	IH6650TD	90 x 45cm	60 x 45cm	50
VEB958	IH6651TD	90 x 90cm	60 x 90cm	24



IH04.05.2017.1

Product Support & Customer Services

Tel: +44 (0) 1453 874130

Email: serviceuk@fannin.eu

www.fannin.eu